Premium Talking Automatic Digital Blood Pressure Monitor

Item # 04-695-001

Please read this guidebook
completely before operating this unit.

English • Spanish

Limited Lifetime Warranty

The warrantor guarantees that its Digital Blood Pressure Monitor will be free from manufacturing defects under normal use for the lifetime from the original purchase date. This warranty covers only normal use, and does not apply to use in any clinical or commercial applications. This warranty does not cover batteries or other power sources that may be provided with or used with the Digital Blood Pressure Monitor. This warranty is voided if the Digital Blood Pressure Monitor product is misused or abused in any manner.

If the Digital Blood Pressure Monitor fails to operate during the lifetime, please call our Customer Care Help Line at 800-622-4714 and we will try to resolve the issue as soon as possible. If the warrantor determines the unit failed to operate due to a manufacturing defect, the unit will be replaced at the option of the warrantor. Replacement is the sole remedy under this limited warranty. This warranty gives you specific legal rights, which vary from state to state. As a condition of this warranty, the enclosed warranty registration card must be completed and sent to us within 10 days of the purchase date or can be filled out online at www.mabisdmi.com.

This warranty constitutes the warrantor's only responsibility and obligation to replace materials or components. We make no other express or implied warranties, arising by operation of law or otherwise, or any warranty of merchantability or fitness for a particular use or purpose whether or not the use or purpose has been disclosed to the warrantor in specifications, drawings or otherwise, and whether or not the warrantor's products are specifically designed and/or manufactured by the warrantor for the buyer's use or purposes, except for the limited warranty stated above. The warrantor will not be responsible for any indirect, incidental, special, consequential, or punitive damages or other loss, including, but not limited to, damage to or loss of other property or equipment and personal injuries, whether to purchaser or others. The warrantor shall in no event be liable to the purchaser for any amount in excess of the cost of repair and/or replacement of the unit.

1931 Norman Drive South • Waukegan, IL 60085
Toll-Free Customer Care Help Line: 1-800-622-4714
Monday - Friday 8:00 am - 4:30 pm CST

Index

Introduction	3
Product Identification & Caution	4
General Blood Pressure Information	5-7
Important Information Before Use	7
World Health Organization (WHO) Indicator	8
Irregular Heartbeat Detection	9
Battery Installation/Replacement	10
Setting the Date/Time	11
Voice Language Setting	12
Applying Your Blood Pressure Cuff	13
Taking Your Blood Pressure Reading	14-15
Recalling Measurements in Memory	15
Description of Display Symbols	16
Troubleshooting	16
Care and Maintenance	17
Product Specifications	17
Blood Pressure Logs	18

Introduction


To achieve the maximum benefit from your blood pressure monitor, we recommend that you first consult with your physician or trained healthcare professional.

Thank you for purchasing an Automatic Digital Blood Pressure Monitor. With proper care and use, your monitor will provide you with many years of reliable readings.

The method of measurement that your Automatic Monitor uses is called the oscillometric method. The monitor detects your blood's movement through the artery in your arm and converts the movements into a digital reading. The oscillometric method does not require a stethoscope, making the monitor easy to use.

Blood pressure readings determined with this device are equivalent to measurements obtained by a trained observer using the cuff/stethoscope auscultation method, within the limits prescribed by the American National Standard for Electronic or Automated Sphygmomanometers.

Product Identification

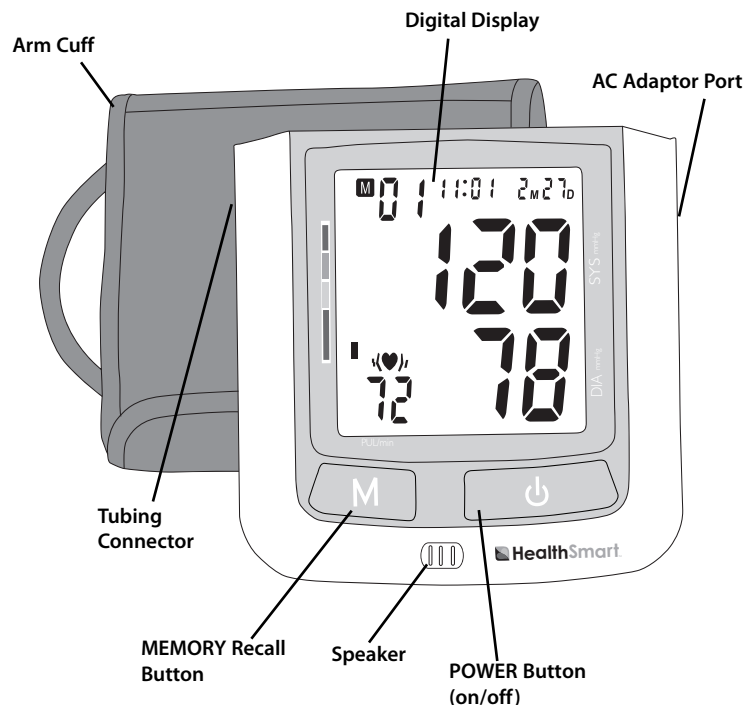CAUTION:

You can stop the inflation or deflation process anytime by pressing the Release button.

General Blood Pressure Information

What is Blood Pressure?

Blood pressure is the pressure that is exerted by blood flowing against the walls of the blood vessels throughout your body.

Your heart, which is the center of the circulatory system, provides the force for the blood to flow or circulate. When your heart contracts, or “beats”, the blood is forced through the blood vessels increasing the pressure. This is the highest pressure in the cycle or what is referred to as **SYSTOLIC** blood pressure. In between beats, your heart relaxes and your blood pressure decreases. This is referred to as **DIASTOLIC** blood pressure.

This complete series of events, which occurs in a single heartbeat, is known as the **CARDIAC CYCLE**.

Your Automatic Blood Pressure Monitor will automatically read your blood pressure and display both systolic and diastolic readings on the screen. Your systolic will be positioned as the upper number and the diastolic reading will be the lower number.

(systolic) **120** / **80** (diastolic)

Blood pressure is measured in millimeters (mm) of mercury (Hg) and is generally recorded with the systolic pressure (120) listed first and the diastolic pressure (80) listed second. The numbers are typically separated by a slash mark (/) as shown.

*Both pressure readings, the **SYSTOLIC** and **DIASTOLIC**, are necessary for a physician to evaluate the status of a patient's blood pressure.*

Please contact your physician for specific information regarding your own blood pressure.

What Influences Blood Pressure?

Many factors such as genetics, age, sex, altitude, physical activity, anxiety, muscular development, certain medications or even the time of day can influence blood pressure. Influences such as sleep or relaxation decrease blood pressure, while anxiety or exercise increase blood pressure.

Why Monitor Your Blood Pressure at Home?

A visit to a physician's office can be a stressful situation for a patient. Anxiety is a known factor in raising blood pressure. This temporary occurrence of elevated blood pressure at the physician's office is commonly referred to as "white coat syndrome".

Whether or not you experience "white coat syndrome", home blood pressure monitoring provides you with the opportunity to supplement your physician's office measurements. These home readings, when taken over a period of time, can show an accurate indication of change. Furthermore, your records can assist your physician in evaluating your health and in making important decisions in the diagnosis and treatment of your condition. Because of this, it is important to take consistent, daily measurements of your blood pressure.

The variations in your individual readings should only be interpreted by your physician or a trained health care professional.

WHO Blood Pressure Classifications


Standards for assessment of high or low blood pressure have been established by the World Health Organization (WHO) as shown on the following chart:

CATEGORY	SYSTOLIC (mmHg)	DIASTOLIC (mmHg)	COLOR INDICATOR
Optimal	<120	<80	GREEN
Normal	120-129	80-84	GREEN
High - Normal	130-139	85-89	GREEN
Stage 1 - Mild Borderline Hypertension	140-159	90-99	YELLOW
Stage 2 - Moderate Hypertension	160-179	100-109	ORANGE
Stage 3 - Severe Hypertension	≥180	≥110	RED
Isolated Systolic Hypertension	≥140	<90	YELLOW
	140-159	<90	ORANGE
	160-179	<90	RED
	≥180	<90	RED

This chart is only a general guideline. Contact your physician or trained healthcare professional to determine your NORMAL blood pressure.

Variations in Blood Pressure

Blood pressure is influenced by many factors and can change from moment to moment. Normally, blood pressure is lowest during sleeping periods and rises during the day. The graph below represents variations in blood pressure shown over a day with measurements taken every 5 minutes.

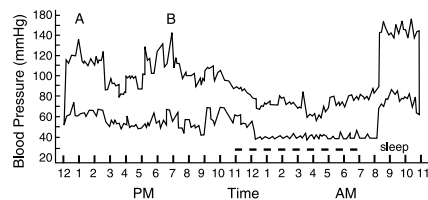
The dotted line represents the sleep period. The rise in blood pressure at 1 p.m. (A in the graph) corresponds to a stressful occurrence and at 7 p.m. (B in the graph) a period of exercise.

Important Information Before Use

1. Blood pressure readings should be interpreted by a physician or trained healthcare professional, who is familiar with your medical history.
2. Perform your measurement in a quiet place. You should be seated and relaxed.
3. Avoid smoking, eating, taking medication, alcohol consumption or physical activity 30 minutes prior to taking a reading. If you are exhibiting signs of stress, avoid taking your measurement until the feeling subsides.
4. Rest 15 minutes prior to taking a reading.
5. Remove any constrictive clothing or jewelry that may interfere with the cuff placement.
6. Keep the monitor stable during measurements. Remain still; do not talk during measurements.
7. Record your daily readings on the chart provided.
8. Take your readings at the same time, each day (or as recommended by your physician).
9. Wait a minimum of 15 minutes between readings. The wait time may vary depending on individual physiological characteristics.
10. This device is intended for adult use.
11. The inflation or deflation process can be stopped by pressing the Release button.

World Health Organization (WHO) Indicator

The WHO Indicator, located on the left side of the display, will identify your reading within a range in accordance with the guidelines set by the World Health Organization (WHO). See the chart below for reference.
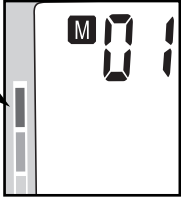
<div><div>SYS</div><div>DIA</div></div>	<div><div>Stage 3 Hypertension</div><div>Systolic pressure: ≥ 180 mmHg -or- Diastolic pressure: ≥ 110 mmHg</div></div>
<div><div>SYS</div><div>DIA</div></div>	<div><div>Stage 2 Hypertension</div><div>Systolic pressure: 160-179 mmHg -or- Diastolic pressure: 100-109 mmHg</div></div>
<div><div>SYS</div><div>DIA</div></div>	<div><div>Stage 1 Hypertension</div><div>Systolic pressure: 140-159 mmHg -or- Diastolic pressure: 90-99 mmHg</div></div>
<div><div>SYS</div><div>DIA</div></div>	<div><div>High Normal</div><div>Systolic pressure: 130-139 mmHg -or- Diastolic pressure: 85-89 mmHg</div></div>
<div><div>SYS</div><div>DIA</div></div>	<div><div>Normal</div><div>Systolic pressure: 120-129 mmHg -or- Diastolic pressure: 80-84 mmHg</div></div>
<div><div>SYS</div><div>DIA</div></div>	<div><div>Optimal</div><div>Systolic pressure: <120 mmHg -and- Diastolic pressure: <80 mmHg</div></div>

Irregular Heartbeat Detection

Your digital blood pressure monitor features Irregular Heartbeat Detection. This feature allows users to accurately monitor blood pressure even if an irregular heartbeat should occur.

When an irregular heartbeat is detected, the IHB icon will appear on the display, Fig. 1.

NOTE: Please consult with your physician or trained healthcare professional for further information regarding an irregular heartbeat if this symbol appears frequently at the end of a reading.

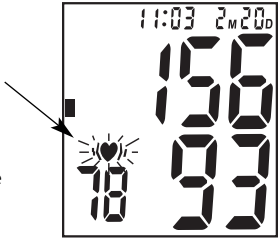
Fig. 1

Battery Installation/Replacement

1. Press down on tab and lift the cover.
2. Insert or replace 4 "AA" batteries into the compartment, matching the indicated polarity symbols, Fig. 1.
3. To close, match up and connect the bottom hooks of the cover, then push the top end.

NOTE: Replace the batteries when the Low Battery symbol appears on the display or when the display remains blank when the POWER button is pressed.


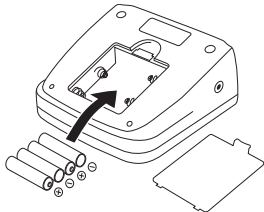
Fig. 1

4. It is recommended to remove the batteries if the unit will not be used for an extended period of time.

AC Power Adapter (optional)

Use the AC adapter for battery-free operation. The AC adapter plugs into the side of the unit as shown, Fig. 2. Please remember to practice electrical safety precautions while using the adapter.

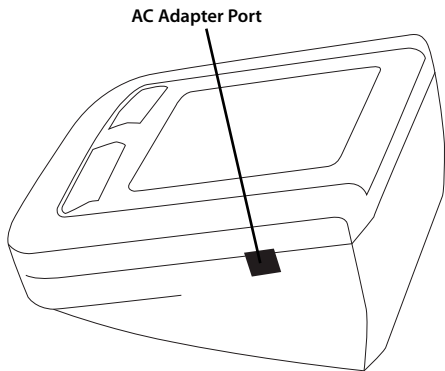Setting the Date/Time

NOTE: When the unit is off, the date and time will still appear at the top of the screen.

1. When the unit is off, press and release both the POWER and MEMORY button at the same time. The month will begin to flash, Fig. 1.
2. Press the MEMORY button to change until the desired number appears.

NOTE: To scroll through the numbers quicker, hold down the MEMORY button.

3. Press the POWER button to set the entry.
4. The date, Fig. 2, followed by the hour and minutes will blink. Repeat steps 2 and 3 to adjust the date/time setting.
5. After setting the minutes, the unit will automatically exit out of the date/time setting mode, Fig. 3.

NOTE: The time is expressed in 24hr military time. The chart below explains the time conversion.

Regular Time (12 hours)	Military Time (24 hours)	Regular Time (12 hours)	Military Time (24 hours)
12 am		12 pm	
Midnight	0	Noon	12
1 am	1	1 pm	13
2 am	2	2 pm	14
3 am	3	3 pm	15
4 am	4	4 pm	16
5 am	5	5 pm	17
6 am	6	6 pm	18
7 am	7	7 pm	19
8 am	8	8 pm	20
9 am	9	9 pm	21
10 am	10	10 pm	22
11 am	11	11 pm	23

Fig. 1

Fig. 2

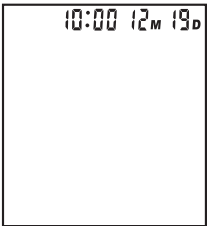
Fig. 3

Voice Language Setting

To set the desired voice language:

1. When the unit is off, you can select the desired voice language by pressing and holding the MEMORY button.
2. Now the screen will cycle through L0, L1 and L2 as shown in Fig. 1.
L0 - voice function is off
L1 - voice is in English
L2 - voice is in Spanish
The monitor will say the corresponding language with the displayed language code.
3. Once the desired language has been selected, release the MEMORY button.
4. After the desired language has been selected, you can adjust the volume by pressing and holding down the MEMORY button until the bars on the side of the screen begin to increase Fig. 2.
5. Continue to hold the MEMORY button down until the volume reaches the desired level, then release the MEMORY button to set.
6. The unit will automatically turn off after the desired volume level has been set.

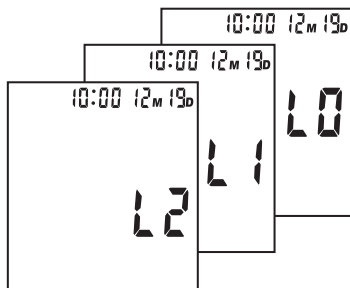
Fig. 1

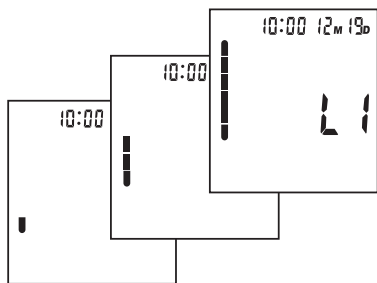
Fig. 2

Applying Your Blood Pressure Cuff

Before applying your blood pressure cuff, be sure you have selected the appropriate cuff:

Standard Cuff fits arm circumference 11-3/4" – 16-1/2" (30-42 cm)

Large Cuff fits arm circumference 16-1/2" – 18-7/8" (42-48 cm)

NOTE: If for any reason you are unable to or should not use your left arm, please modify the following instructions and apply the cuff to your right arm. Your physician can tell you which arm is best for you to use.

1. Remove any constrictive clothing or jewelry that may interfere with the cuff placement.
2. Be seated with your feet flat on the floor.
3. Position the cuff on a solid surface with the tubing facing UP and AWAY from you, Fig. 1. The metal bar on the cuff should be to the left of the tubing.
4. Widen or open the cuff by pulling or rolling the bottom of the cuff towards the right, Fig. 2. This should open the cuff, creating a cylinder. Do not extend the cuff beyond the metal bar.

Fig. 1

Fig. 2

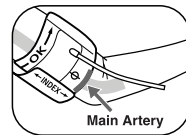
Fig. 3

5. Insert your arm into the cuff (cylinder). Position the (Φ) mark over the main artery (on the inside of your arm), Fig. 3.
6. The bottom edge of the cuff should be positioned approximately one inch ABOVE the elbow joint, Fig. 4.
7. Reaching underneath your left arm with your right hand, pull the end of the cuff towards your body to tighten the cuff, Fig. 5. Wrap and secure the cuff making sure that the (Φ) mark remains as shown, Fig. 3.
8. The cuff should fit comfortably, yet snugly around your arm. You should be able to insert two fingers between your arm and the cuff.

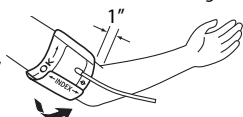
Fig. 4


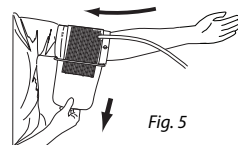
Fig. 5

NOTE: The cuff size is suitable for use when the vertical "Index" mark is within the horizontal "OK" range. A different sized cuff is needed if the "Index" mark is outside of the "OK" range.

Taking Your Blood Pressure Reading

Proceed only after reading the previous sections of this manual.


1. Position the monitor on a flat, stable surface with the digital display panel in view. Insert the cuff tubing connector into the left side of the monitor, Fig. 1.

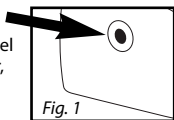2. Rest your elbow on a solid surface with your palm facing upward. Elevate your arm so that the cuff is at the same level as your heart, Fig. 2. Relax your left hand.

3. Press the POWER button to turn the unit on. Select the memory bank (1 or 2) to store the reading in by pressing the MEMORY button. Press POWER to confirm your selection.

4. The monitor will begin to guide you with voice prompts if the voice function is turned on. This unit will run a self-test, Fig. 3.

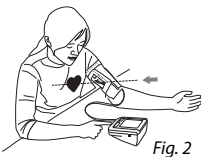5. This unit will show the last measured reading, Fig. 4, or if there is no previous measurement, "0" as shown in Fig. 5.

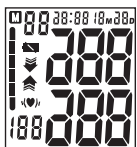
Fig. 3

Fig. 4

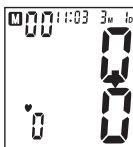
Fig. 5

6. When the display appears as shown, Fig. 6, the monitor is ready and will begin to guide you with voice prompts if the voice function is on. The cuff will automatically begin to inflate.

7. Once the cuff inflates to the appropriate pressure level, the measurement will begin.

8. The unit will begin to automatically deflate and the numbers appearing on the display will decrease.

9. When all the air is deflated, the measurement is complete. The unit will display your blood pressure (systolic and diastolic) measurements and pulse rate, Fig. 7. With the voice function on (L1 or L2), the monitor will say the measurement with the language you selected according to the classification standard of WHO (World Health Organization). Record your readings in a chart.

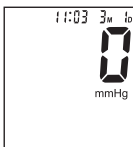
Fig. 6

Fig. 7

10. Your measurements will automatically be saved in memory.

11. Press the POWER button to turn the power off, or the unit will automatically shut off after approximately 1 minute.

12. Disconnect the tubing from the monitor prior to storing.

NOTE: This monitor inflates the cuff to approximately 190 mmHg or 40 mmHg higher than the last systolic pressure stored in memory. If the system detects that the cuff needs additional pressure, it will automatically re-inflate to approximately 40 mmHg higher than the last pressure level.

Recalling Measurements in Memory

You can recall up to 120 total measurements in memory (60 per memory bank). These measurements can be shared with your physician or trained healthcare professional.

1. Press and release the MEMORY button. The unit will display either U1 or U2 depending upon the user bank selected. It will also show how many measurements are stored within that bank.

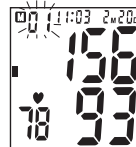
Fig. 1

2. To select the appropriate user bank press the POWER button to switch between U1 and U2.

3. Once the desired user memory bank is selected, press the MEMORY button again to scroll through the saved measurements, Fig. 1. After you've scrolled through all the readings, the unit will automatically turn off.

NOTE: Each memory bank can store up to 60 readings. When the number of readings exceeds 60, the oldest data will be replaced with the new record.

To erase all stored measurements:

1. Press and release the MEMORY button. The unit will display either U1 or U2 depending upon the user bank selected and how many measurements are stored within that bank.
2. Select the desired user bank by pressing and releasing the MEMORY button.
3. Press and hold down the Memory button until two horizontal lines flash, Fig. 2. This indicates that all measurements within the selected memory bank have been erased.

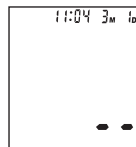
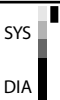
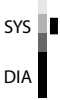
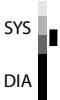
Fig. 2

Description of Display Symbols

	Indicates a measurement within Stage 3 Hypertension.		Shows the pulse rate per minute.
			Occurs when an error was made during measurement.
	Indicates a measurement within Stage 2 Hypertension.		Appears when there is an irregular heartbeat detected during measurement.
			Appears when the cuff is unstable or it has too much air left in it. Once the arrow stops flashing, the unit is ready to take a measurement.
	Indicates a measurement within Stage 1 Hypertension.		Appears when the unit is ready or a measurement error has occurred. Restart the measurement if an error occurs.
			Appears when recalling a measurement from memory.
	Appears when batteries should be replaced.		

Troubleshooting

If any abnormality occurs during use, please check and correct the following:

Condition	Correction
Display is blank when power is on.	Check and correct the polarity of the installed batteries. Reinstall or replace batteries.
Measurement incomplete or abnormally low or high values displayed.	Review and follow "Applying Your Blood Pressure Cuff" and "Taking Your Blood Pressure Reading" sections. Replace the batteries.
EE symbol displays.	Do not talk or move during the measurement. Review "Applying Your Blood Pressure Cuff" section.
Measurements are different from those typically measured by physician or every measurement reading is different.	Remember, blood pressure readings are influenced by physical and mental conditions and/or even the time of day. Daily measurements should be interpreted by your physician.

Care and Maintenance

1. Disconnect the cuff tubing from the monitor and remove batteries prior to storing.
2. Only use a soft, dry cloth to clean your blood pressure monitor.
3. Avoid using any types of liquids on the monitor or cuff.
4. Do not store the unit where it will be exposed to direct sunlight, dust or humidity.
5. Avoid extreme temperatures.
6. Never disassemble the monitor or cuff.
7. Dropping or subjecting your blood pressure monitor to strong shocks should be avoided.

Specifications

Name and model number	Talking Automatic Digital Blood Pressure Arm Monitor, 04-695-001
Display system	Digital display/LCD
Measuring method	Oscillometric
Power source	4 "AA" batteries
Measuring range	Pressure: 0-300 mmHg Pulse: 40-180 beats/minute
Accuracy	Pressure: ± 3 mmHg Pulse: $\pm 5\%$ of reading
Inflation	Automatic Inflation
Deflation	Automatic Pressure Release Valve/Manual fast release
Memory	Built-in memory enabling display of up to 120 measurements
Automatic shut-off	Approximately 1 minute after last button operation
Operation environment	Temperature 41°F – 104°F (5°C – 40°C) Humidity < 90%
Storage environment	Temperature -4°F – 131°F (-20°C – 55°C) Humidity < 95%
Monitor dimensions	5" x 5-1/8" x 2-3/8"
Weight	11.11 oz. (without batteries)
Arm circumference range	Adult cuff fits arm circumference: 11-3/4" – 16-1/2" (30-42 cm) Large cuff fits arm circumference 16-1/2" – 18-7/8" (42-48 cm)
Accessories	Cuff, batteries, case, detailed guidebook, quick start guide, blood pressure log

Specifications are subject to change without notice

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