

JOLLY DIGITAL BLOOD PRESSURE MONITOR

Operation Guide




ATTENTION:
The operators must carefully read and completely understand the present manual before using the product.

SPECIFICATIONS

- Product name: Arm Blood Pressure Monitor
- Model: 32777 / KD-5920
- Classification: Internally powered, Type BF applied part, IP20, No AP or APG, Continuous operation
- Machine size: Approx.150 mm × 95 mm × 41 mm (5 29/32" × 3 3/4" × 1 5/8")
- Cuff circumference: 22-30 cm (8-11"), 30-42 cm (11-16") optional, 42-48 cm (16-18.2") optional
- Weight: Approx. 235 g (8 9/32 oz.) (exclude batteries)
- Measuring method: Oscillometric method, automatic inflation and measurement
- Memory volume: 4×30 times with time and date stamp
- Power source: batteries: 4×1.5 V  SIZE AAA
- Measurement range:
 - Cuff pressure: 0-300 mmHg
 - Systolic: 60-260 mmHg
 - Diastolic: 40-199 mmHg
 - Pulse rate: 40-180 beats/minute
- Accuracy:
 - Pressure: ±3 mmHg
 - Pulse rate: ±5%
- Environmental temperature for operation: 10°C~40°C (50°F~104°F)
- Environmental humidity for operation: ≤85% RH

2. Clock and date adjustment

- Once you install the battery or turn off the monitor, it will enter Clock Mode, and LCD will display time and date.
- While the monitor is in Clock Mode, press the "START" and "MEM" button simultaneously for two seconds, the month will blink at first. Press the button "START" repeatedly, the day, hour and minute will blink in turn. While the number is blinking, press the button "MEM" to increase the number. Keep on pressing the button "MEM", the number will increase fast.

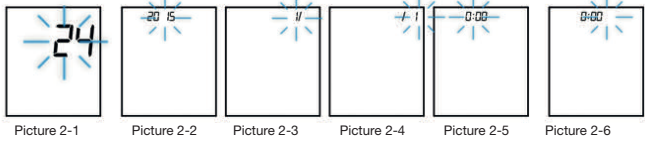

- You can turn off the monitor by pressing "START" button when the minute is blinking, then the time and date is confirmed.
- The monitor will turn off automatically after 1 minute of no operation with the time and date unchanged.
- Once you change the batteries, you should readjust the time and date.


IMPORTANT INFORMATION	3
CONTENTS AND DISPLAY INDICATORS	4
INTENDED USE	5
CONTRAINDICATION	5
PRODUCT DESCRIPTION	6
SPECIFICATIONS	7
NOTICE	9
SETUP AND OPERATING PROCEDURES	12
1. Battery loading	12
2. Clock and date adjustment	13
3. Connecting the cuff to the monitor	14
4. Applying the cuff	14
5. Body posture during measurement	15
6. Taking your blood pressure reading	16
7. Displaying stored results	18
8. Deleting measurements from the memory	20
9. Assessing high blood pressure for adults	20
10. Technical alarm description	22
11. Troubleshooting (1)	23
12. Troubleshooting (2)	24
MAINTENANCE	26
EXPLANATION OF SYMBOLS ON UNIT	28
ELECTROMAGNETIC COMPATIBILITY INFORMATION	30
WARRANTY INFORMATION	35

- Environmental temperature for storage and transport: -20°C~50°C (-4°F~122°F)
- Environmental humidity for storage and transport: ≤85%RH
- Environmental pressure: 80kPa-105kPa
- Battery life: Approx 100 times.
- All components belonging to the pressure measuring system, including: Pump,Valve, LCD, Cuff, Sensor

Note: These specifications are subject to change without notice.

3. Connecting the cuff to the monitor



Insert the cuff tubing connector into the socket in the left side of the monitor. Make certain that the connector is completely inserted to avoid air leakage during blood pressure measurements.

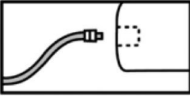

 Avoid compression or restriction of the connection tubing during measurement, which may cause inflation error, or harmful injury due to continuous cuff pressure.

4. Applying the cuff

- Pulling the cuff end through the medial loop (the cuff is packaged like this already), turn it outward (away from your body) and tighten it and close the Velcro fastener.
- Place the cuff around a bare arm 1-2cm above the elbow joint.
- While seated, place palm upside in front of you on a flat surface such as a desk or table. Position the air tube in the middle of your arm in line with your middle finger.

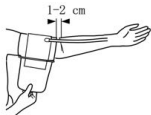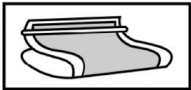

IMPORTANT INFORMATION

Normal blood pressure fluctuation

All physical activity, excitement, stress, eating, drinking, smoking, body posture and many other activities or factors (including taking a blood pressure measurement) will influence blood pressure value. Because of this, it is mostly unusual to obtain identical multiple blood pressure readings.

Blood pressure fluctuates continually ----- day and night. The highest value usually appears in the daytime and lowest one usually at midnight. Typically, the value begins to increase at around 3:00AM, and reaches to highest level in the daytime while most people are awake and active.

Considering the above information, it is recommended that you measure your blood pressure at approximately the same time each day.

Too frequent measurements may cause injury due to blood flow interference, please always relax a minimum of 1 to 1.5 minutes between measurements to allow the blood circulation in your arm to recover. It is rare that you obtain identical blood pressure readings each time.

NOTICE

- Read all of the information in the operation guide and any other literature in the box before operating the unit.
- Stay still, calm and rest for 5 minutes before blood pressure measurement.
- The cuff should be placed at the same level as your heart.
- During measurement, neither speak nor move your body and arm.
- Measuring on same arm for each measurement.
- Please always relax at least 1 or 1.5 minutes between measurements to allow the blood circulation in your arm to recover. Prolonged over-inflation (cuff pressure exceed 300 mmHg or maintained above15 mmHg for longer than 3 minutes) of the bladder may cause ecchymoma of your arm.
- Consult your physician if you have any doubt about below cases:
 - The application of the cuff over a wound or inflammation diseases;
 - The application of the cuff on any limb where intravascular access or therapy, or an arterio-venous (A-V) shunt, is present;
 - The application of the cuff on the arm on the side of a mastectomy;
 - Simultaneously used with other monitoring medical equipments on the same limb;
 - Need to check the blood circulation of the user.
-  This Electronic Sphygmomanometers is designed for adults and should never be used on infants or young children. Consult your physician or other health care professionals before use on older children.

- The cuff should fit comfortably, yet snugly around your arm. You should be able to insert one finger between your arm and the cuff.

Note:

- Please refer to the cuff circumference range in "SPECIFICATIONS" to make sure that the appropriate cuff is used.
- Measuring on same arm each time.
- Do not move your arm, body, or the monitor and do not move the rubber tube during measurement.
- Stay quiet, calm for 5 minutes before blood pressure measurement.
- Please keep the cuff clean. If the cuff becomes dirty, remove it from the monitor and clear it by hand in a mild detergent, then rinse it thoroughly in cold water. Never dry the cuff in clothes dryer or iron it. Clean the cuff after the usage of every 200 times is recommended.


5. Body posture during measurement

Sitting Comfortably Measurement

- Be seated with your feet flat on the floor, and don't cross your legs.
- Place palm upside in front of you on a flat surface such as a desk or table.
- The middle of the cuff should be at the level of the right atrium of the heart.



CONTENTS AND DISPLAY INDICATORS

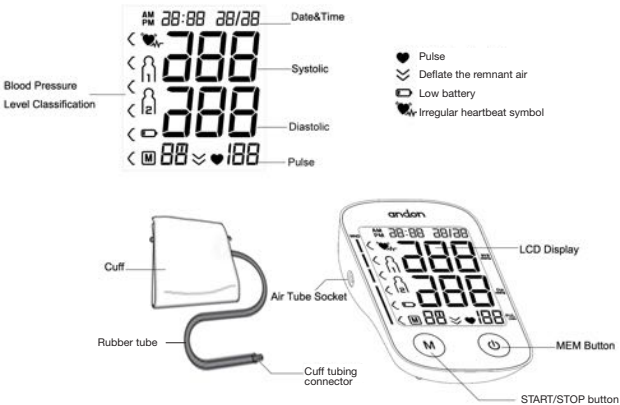

- Do not use this unit in a moving vehicle, This may result in erroneous measurement.
- Blood pressure measurements determined by this monitor are equivalent to those obtained by a trained observer using the cuff/stethoscope auscultation method, within the limits prescribed by the American National Standard Institute, Electronic or automated sphygmomanometers.
- Information regarding potential electromagnetic or other interference between the blood pressure monitor and other devices together with advice regarding avoidance of such interference please see part ELECTROMAGNETIC COMPATIBILITY INFORMATION.
- If Irregular Heartbeat (IHB) brought by common arrhythmias is detected in the procedure of blood pressure measurement, a signal of  will be displayed. Under this condition, the Electronic Sphygmomanometers can keep function, but the results may not be accurate, it's suggested that you consult with your physician for accurate assessment.
There are 2 conditions under which the signal of IHB will be displayed:
 - The coefficient of variation (CV) of pulse period >25%.
 - The difference of adjacent pulse period≥0.14s, and the number of such pulse takes more than 53 percentage of the total number of pulse.
- Please do not use the cuff other than supplied by the manufacturer, otherwise it may bring biocompatible hazard and might result in measurement error.


Lying Down Measurement

- Lie on your back.
- Place your arm straight along your side with your palm upside.
- The cuff should be placed at the same level as your heart.



6. Taking your blood pressure reading

- After applying the cuff and your body is in a comfortable position, press the "START" button. A beep is heard and all display characters are shown for self-test. See Picture 7-1. Please contact the service center if segment is missing.

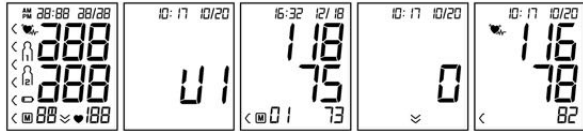

Picture 7-1 Picture 7-2 Picture 7-3 Picture 7-4 Picture 7-5

INTENDED USE

Fully Automatic Electronic Sphygmomanometer is for use by medical professionals or at home and is a non-invasive blood pressure measurement system intended to measure the diastolic and systolic blood pressures and pulse rate of an adult individual by using a non-invasive technique in which an inflatable cuff is wrapped around the upper arm. The cuff circumference is limited to 22cm-42cm (approx. 8 21/32"-16 17/32").

CONTRAINDICATION

 It is inappropriate for people with serious arrhythmia to use this Electronic Sphygmomanometer.

-  The monitor might not meet its performance specifications or cause safety hazard if stored or used outside the specified temperature and humidity ranges in specifications.
-  Please do not share the cuff with other infective person to avoid cross-infection.
- This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:
 - Reorient or relocate the receiving antenna.
 - Increase the separation between the equipment and receiver.
 - Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
 - Consult the dealer or an experienced radio/TV technician for help.

- Then the current memory bank (U1, U2, U3 or U4) is shown. See picture 7-2. Press "MEM" button to change over to other bank. Confirm your selection by pressing "START" button. The current bank can also be confirmed automatically after 5 seconds with no operation.
- If the monitor has stored results, the LCD will momentarily display the most recent one. If no result has been stored, zero will appear on LCD. See Picture 7-3.
- Then the monitor starts to seek zero pressure. See picture 7-4.
- Then the monitor inflates the cuff until sufficient pressure has built up for a measurement. Then the monitor slowly releases air from the cuff and carries out the measurement. Finally the blood pressure and pulse rate will be calculated and displayed on the LCD screen. The blood pressure classification indicator and Irregular heartbeat symbol (if any) will blink on the screen. If the voice function is on, it will announce the measurement result. The result will be automatically stored in the monitor. See Picture 7-5.
- After measurement, the monitor will turn off automatically after 1 minute of no operation. Alternatively, you can press the "START" button to turn off the monitor manually.
- During measurement, you can press the "START" button to turn off the monitor manually.

Note: Please consult a health care professional for interpretation of pressure measurements.

PRODUCT DESCRIPTION

Based on Oscillometric methodology and silicon integrated pressure sensor, blood pressure and pulse rate can be measured automatically and non-invasively. The LCD display will show blood pressure and pulse rate. The most recent 4x30 measurements can be stored in the memory with date and time stamp. The voice function will ease the operation. The Electronic Sphygmomanometers corresponds to the below standards: IEC 60601-1:2005 +A1:2012(E)/EN 60601-1:2006/A11: 2011 (Medical electrical equipment -- Part 1: General requirements for basic safety and essential performance), IEC60601-1-2:2007/EN 60601-1-2:2007 /AC:2010 (Medical electrical equipment -- Part 1-2: General requirements for basic safety and essential performance - Collateral standard: Electromagnetic compatibility - Requirements and tests), IEC80601-2-30:2009+AMD1:2013/EN 60601-2-30:2010/A1:2015(Medical electrical equipment -- Part 2-30: Particular requirements for the basic safety and essential performance of automated non-invasive sphygmomanometers)EN 1060-1: 1995 + A2: 2009 (Non-invasive sphygmomanometers - Part 1: General requirements), EN 1060-3: 1997 + A2: 2009 (Non-invasive sphygmomanometers - Part 3: Supplementary requirements for electro-mechanical blood pressure measuring systems).

SETUP AND OPERATING PROCEDURES

1. Battery loading

- Open battery cover at the back of the monitor.
- Load four "AAA" size batteries. Please pay attention to polarity.
- Close the battery cover.


When LCD shows battery symbol  , replace all batteries with new ones. Rechargeable batteries are not suitable for this monitor. Remove the batteries if the monitor will not be used for a month or more to avoid relevant damage of battery leakage.

 Avoid the battery fluid to get in your eyes. If it should get in your eyes, immediately rinse with plenty of clean water and contact a physician.

 The monitor, the batteries and the cuff, must be disposed of according to local regulations at the end of their usage.


7. Displaying stored results

- After the measurement, you can review the measurements in the current memory bank by pressing button "MEM". Now the LCD displays the amount of the results in the current bank.

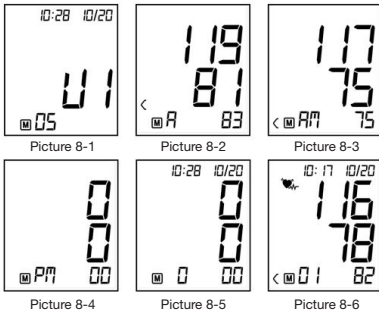


Picture 8-1 Picture 8-2 Picture 8-3 Picture 8-4 Picture 8-5 Picture 8-6

- b. Alternatively, press “MEM” button in Clock Mode to display the stored results. The current memory bank will blink and the amount of results in this bank will be displayed. Press “START” button to change over to other bank. Confirm your selection by pressing “MEM” button. The current bank can also be confirmed automatically after 5 seconds with no operation. See picture 8-1.
- c. After selecting the memory bank, the LCD will display the average values in this bank. See picture 8-2&8-3&8-4. If no result stored, LCD will show zeros. See picture 8-5. Press “MEM” button to switch to next scene.
- d. When the average is displayed, press the “MEM” button, the most recent result will be displayed. See picture 8-6. Followed by, the blood pressure and pulse rate will be shown separately. Irregular heartbeat symbol (if any) will blink. Press “MEM” button again to review the next result. In this way, repeatedly pressing the “MEM” button displays the respective results measured previously.
- e. When displaying the stored results, the monitor will turn off automatically after 1 minute of no operation. You can also press the button “START” to turn off the monitor manually.

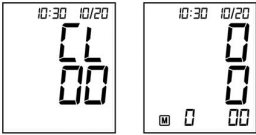

9. Assessing high blood pressure for adults

The following guidelines for assessing high blood pressure (without regard to age or gender) have been established by the World Health Organization (WHO). Please note that other factors (e.g. diabetes, obesity, smoking, etc.) need to be taken into consideration. Consult with your physician for accurate assessment, and never change your treatment by yourself.

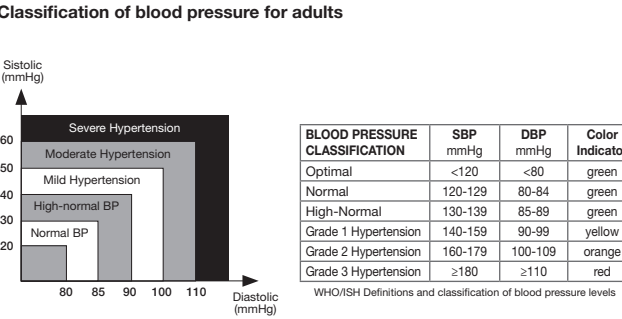

Note: It is not intended to provide a basis of any type of rush toward emergency conditions/diagnosis based on the color scheme and that the color scheme is meant only to discriminate between the different levels of blood pressure.

10. Technical alarm description

The monitor will show ‘Hi’ or ‘Lo’ as technical alarm on LCD with no delay if the determined blood pressure (systolic or diastolic) is outside the rated range specified in part SPECIFICATIONS. In this case, you should consult a physician or check if your operation violated the instructions.

The technical alarm condition (outside the rated range) is preset in the factory and cannot be adjusted or inactivated. This alarm condition is assigned as low priority according to IEC 60601-1-8.

The technical alarm is non-latching and need no reset. The signal displayed on LCD will disappear automatically after about 8 seconds.

PROBLEM	POSSIBLE CAUSE	SOLUTION
LCD Display shows abnormal result	The cuff position was not correct or it was not properly tightened	Apply the cuff correctly and try again
	Body posture was not correct during testing	Review the “BODY POSTURE DURING MEASUREMENT” sections of the instructions and re-test
	Speaking, arm or body movement, angry, excited or nervous during testing	Re-test when calm and without speaking or moving during the test
	Irregular heartbeat (arrhythmia)	It is inappropriate for people with serious arrhythmia to use this Electronic Sphygmomanometer

PROBLEM	POSSIBLE CAUSE	SOLUTION
LCD Display shows battery symbol	Low Battery	Change the batteries
LCD shows “Er 0”	Pressure system is unstable before measurement	Don’t move and try again
LCD shows “Er 1”	Fail to detect systolic pressure	
LCD shows “Er 2”	Fail to detect diastolic pressure	
LCD shows “Er 3”	Pneumatic system blocked or cuff is too tight during inflation	Apply the cuff correctly and try again
LCD shows “Er 4”	Pneumatic system leakage or cuff is too loose during inflation	

PROBLEM	POSSIBLE CAUSE	SOLUTION
LCD shows “Er 5”	Cuff pressure above 300mmHg	Measure again after five minutes. If the monitor is still abnormal, please contact the local distributor or the factory
LCD shows “Er 6”	More than 160 seconds with cuff pressure above 15 mmHg	
LCD shows “Er 7”	Inner memory error	
LCD shows “Er 8”	Device parameter checking error	
LCD shows “Er A”	Pressure sensor parameter error	Take out batteries for five minutes, and then reinstall all batteries
No response when you press button or load battery	Incorrect operation or strong electromagnetic interference	

- MAINTENANCE**
1. Do not drop this monitor or subject it to strong impact.
2. Avoid high temperature and solarization. Do not immerse the monitor in water as this will result in damage to the monitor.
3. If this monitor is stored near freezing, allow it to acclimate to room temperature before use.
4. Do not attempt to disassemble this monitor.
5. If you do not use the monitor for a long time, please remove the batteries.
6. It is recommended the performance should be checked every 2 years or after repair. Please contact the service center.
7. Clean the monitor with a dry, soft cloth or a soft cloth squeezed well after moistened with water, diluted disinfectant alcohol, or diluted detergent.
8. No component can be maintained by user in the monitor. The circuit diagrams, component part lists, descriptions, calibration instructions, or other information which will assist the user's appropriately qualified technical personnel to repair those parts of equipment which are designated repairably can be supplied.
9. The monitor can maintain the safety and performance characteristics for a minimum of 10,000 measurements or three years, and the cuff integrity is maintained after 1,000 open–close cycles of the closure.

10. It is recommended the cuff should be disinfected 2 times every week if needed (For example, in hospital or in clinique). Wipe the inner side (the side contacts skin) of the cuff by a soft cloth squeezed after moistened with Ethyl alcohol (75-90%), then dry the cuff by airing.

EXPLANATION OF SYMBOLS ON UNIT

Symbol for “THE OPERATION GUIDE MUST BE READ”
(The sign background colour: blue. The sign graphical symbol: white)

Symbol for “WARNING”

Symbol for “TYPE BF APPLIED PARTS”
(The cuff is type BF applied part)

Symbol for “ENVIRONMENT PROTECTION – Waste electrical products should not be disposed of with household waste. Please recycle where facilities exist. Check with your local Authority or retailer for recycling advice”.

Symbol for “MANUFACTURER”

Symbol for “COMPILES WITH MDD93/42/EEC REQUIREMENTS”

Symbol for “DATE OF MANUFACTURE”

Symbol for “EUROPEAN REPRESENTATION”

- SN Symbol for “SERIAL NUMBER”
- Symbol for “KEEP DRY”
- Symbol for “KEEP AWAY FROM SUNLIGHT”
- Symbol for “PRODUCT CODE”
- Symbol for “LOT NUMBER”

ELECTROMAGNETIC COMPATIBILITY INFORMATION

Guidance and manufacture's declaration - electromagnetic emissions		
The 32777 is intended for use in the electromagnetic environment specified below. The customer or the user of the 32777 should assure that it is used in such an environment.		
Emissions test	Compliance	Electromagnetic environment - guidance
RF emissions CISPR 11	Group 1	The 32777 uses RF energy only for its internal function. Therefore, its RF emissions are very low and are not likely to cause any interference in nearby electronic equipment.
RF emissions CISPR 11	Class B	The 32777 is suitable for use in all establishments other than domestic and those directly connected to the public low-voltage power supply network that supplies buildings used for domestic purposes.
Harmonic emissions IEC 61000-3-2	Not applicable	
Voltage fluctuations/ flicker emissions IEC 61000-3-3	Not applicable	

Guidance and manufacturer's declaration - electromagnetic immunity			
The 32777 is intended for use in the electromagnetic environment specified below. The customer or the user of the 32777 should assure that it is used in such an environment.			
IMMUNITY TEST	IEC 60601test level	Compliance level	Electromagnetic environment - guidance
Electrostatic discharge (ESD) IEC 61000-4-2	± 6 kV contact ± 8 kV air	± 6 kV contact ± 8 kV air	Floors should be wood, concrete or ceramic tile. If floors are covered with synthetic material, the relative humidity should be at least 30%.
Power frequency (50/60 Hz) magnetic field IEC 61000-4-8	3 A/m	3 A/m	Power frequency magnetic fields should be at levels characteristic of a typical location in a typical commercial or hospital environment.

Guidance and manufacturer's declaration - electromagnetic immunity			
The 32777 is intended for use in the electromagnetic environment specified below. The customer or the user of the 32777 should assure that it is used in such an environment.			
IMMUNITY TEST	IEC 60601 test level	Compliance level	Electromagnetic environment - guidance
Radiated RF IEC 61000-4-3	3 V/m 80 MHz to 2.5 GHz	3 V/m	Portable and mobile RF communications equipment should be used no closer to any part of the 32777, including cables, than the recommended separation distance calculated from the equation applicable to the frequency of the transmitter. Recommended separation distance $d = 2.3 \sqrt{P}$ 800 MHz to 2.5 GHz

			Where P is the maximum output power rating of the transmitter in watts (W) according to the transmitter manufacturer and d is the recommended separation distance in meters (m). Field strengths from fixed RF transmitters, as determined by an electromagnetic site survey, ^a should be less than the compliance level in each frequency range. ^b Interference may occur in the vicinity of equipment marked with the following symbol:
NOTE 1 At 80 MHz and 800 MHz, the higher frequency range applies. NOTE 2 These guidelines may not apply in all situations. Electromagnetic propagation is affected by absorption and reflection from structures, objects and people.			
a Field strengths from fixed transmitters, such as base stations for radio (cellular/cordless) telephones and land mobile radios, amateur radio, AM and FM radio broadcast and TV broadcast cannot be predicted theoretically with accuracy. To assess the electromagnetic environment due to fixed RF transmitters, an electromagnetic site survey should be considered. If the measured field strength in the location in which the 32777 is used exceeds the applicable RF compliance level above, the 32777 should be observed to verify normal operation. If abnormal performance is observed, additional measures may be necessary, such as re-orienting or relocating the 32777. b Over the frequency range 150 kHz to 80 MHz, field strengths should be less than 3 V/m.			

Recommended separation distances between portable and mobile RF communications equipment and the 32777			
The 32777 is intended for use in an electromagnetic environment in which radiated RF disturbances are controlled. The customer or the user of the 32777 can help prevent electromagnetic interference by maintaining a minimum distance between portable and mobile RF communications equipment (transmitters) and the 32777 as recommended below, according to the maximum output power of the communications equipment.			
Rated maximum output power of transmitter W	Separation distance according to frequency of transmitter m		
	150 kHz to 80 MHz $d = 1.2 \sqrt{P}$	80 MHz to 800 MHz $d = 1.2 \sqrt{P}$	800 MHz to 2.5 GHz $d = 2.3 \sqrt{P}$
0,01	0,12	0,12	0,23
0,1	0,38	0,38	0,73
1	1,2	1,2	2,3
10	3,8	3,8	7,3
100	12	12	23
For transmitters rated at a maximum output power not listed above, the recommended separation distance d in meters (m) can be determined using the equation applicable to the frequency of the transmitter, where P is the maximum output power rating of the transmitter in watts (W) according to the transmitter manufacturer. NOTE 1 At 80 MHz and 800 MHz, the separation distance for the higher frequency range applies. NOTE 2 These guidelines may not apply in all situations. Electromagnetic propagation is affected by absorption and reflection from structures, objects and people.			

Disposal: The product must not be disposed of along with other domestic waste. The users must dispose of this equipment by bringing it to a specific recycling point for electric and electronic equipment. For further information on recycling points contact the local authorities, the local recycling center or the shop where the product was purchased. If the equipment is not disposed of correctly, fines or penalties may be applied in accordance with the national legislation and regulations.

GIMA WARRANTY CONDITIONS
Congratulations for purchasing a GIMA product. This product meets high qualitative standards both as regards the material and the production.
The warranty is valid for 12 months from the date of supply of GIMA.
During the period of validity of the warranty, GIMA will repair and/or replace free of charge all the defected parts due to production reasons.
Labor costs and personnel traveling expenses and packaging not included.
All components subject to wear are not included in the warranty.
The repair or replacement performed during the warranty period shall not extend the warranty. The warranty is void in the following cases: repairs performed by unauthorized personnel or with non-original spare parts, defects caused by negligence or incorrect use. GIMA cannot be held responsible for malfunctioning on electronic devices or software due to outside agents such as: voltage changes, electro-magnetic fields, radio interferences, etc.
The warranty is void if the above regulations are not observed and if the serial code (if available) has been removed, cancelled or changed.
The defected products must be returned only to the dealer the product was purchased from. Products sent to GIMA will be rejected.

32777 / KD-5920

ANDON HEALTH CO., LTD.
No. 3 JinPing Street, Ya An Road, Nankai District,
Tianjin 300190, China
Made in China

Lotus Global Co., Ltd.
1 Four Seasons Terrace, West Drayton,
Middlesex, London, UB7 9GG, United Kingdom

M-32777-GB-Rev0.1.1.17