GIMA

PROFESSIONAL MEDICAL PRODUCTS

Gima S.p.A. - Via Marconi, 1 - 20060 Gessate (MI) Italy
Italia: tel. 199 400 401 - fax 199 400 403
Export: tel. +39 02 953854209/221/225 fax +39 02 95380056
gima@gimaitaly.com - export@gimaitaly.com
www.gimaitaly.com

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
NAME OF PARTS AND FUNCTIONS

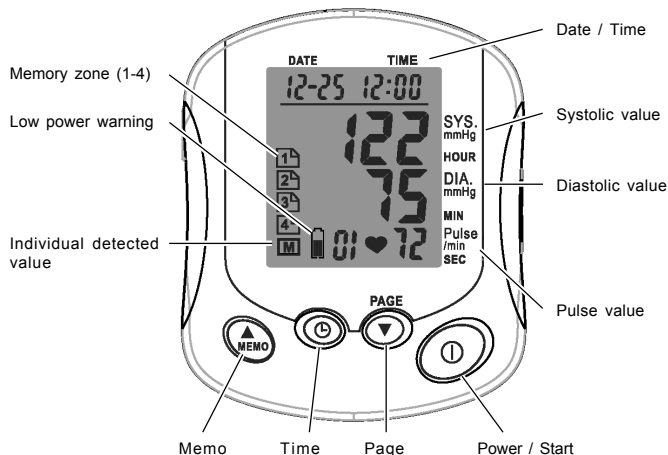Function

- Automatic determines ideal cuff inflation. 4-person memory zones, each 25-reading with measuring date/time.
- Super LCD display easy to read
- Auto shut-off in 3 minutes, power saving
- Operated by 4 x "AA" batteries

POWER/START: power on/off or start measuring/measure again

PAGE: press "PAGE" for choosing GROUP (1~4)

MEMO: Recall different previous values

TIME: for time adjustment, press "▲ TIME" about 5 seconds until the digits of year/month/date/time blinking, start to press "MEMO" for time forward; press "▼ PAGE" for time backward.



HOW TO MEASURE BLOOD PRESSURE

1. To insert the batteries

- Remove the battery cover in accordance with the direction of arrow.
- Insert 4 x "AA" batteries in accordance with the indication of polarities on the batteries.
- Slide the battery cover to close.

2. Measurement

Connect the tube of arm cuff to the unit properly.

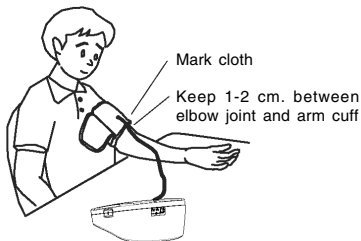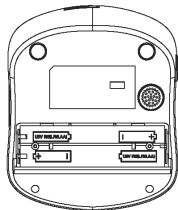1. To apply the arm cuff.

- (1) Turn the palm of your hand upward. Place mark cloth inside your arm and match the end of the arm cuff to approximately 1-2 cm from the flexed side of the elbow joint of your arm (nipple level) and wrap the arm cuff on the arm.
- (2) Hold the end of the Arm Cuff and wrap it up by gently pulling so that there is no gap between the arm cuff and the arm.
- (3) Press "Power" button to measure blood pressure monitor after wrapped arm cuff.

2. To start the measurement

Note: In order to get correct measurement results and to avoid unnecessary erroneous measurement, you should keep your arm hand and fingers be still during measurement, be seated and not talk.

- (A) Push "START" button, the cuff will inflate automatically
- (B) After main unit detecting Systolic/Diastolic and Pulse Values etc. all of the measured values will be displayed on the LCD panel.
- (C) Release Arm Cuff from your arm while the cuff deflates completely.

3. To see the previous measured values

The main unit automatically stores the previous measured values.

(Even after batteries are replaced, the previous values are stored in one minute.)

- (A) Press "MEMO" button, then you can recall previous values for the first person's measured value.
- (B) If you press "MEMO" BUTTON SUCCESSIVELY then you can read next values.
- (C) Press "PAGE" button, then you can read second, third or forth person's measured values.
- (D) Stored capacity is up to 25 sets each page.

4. To finish the measurement

After deflation, the measurement of blood pressure is completed, if you want to measure again, just press "Power" button again, the unit will shut-off and you should press "Power" button again to restart the power.

HOW TO SET TIME AND DATE

Generally, LCD panel displays Month/Date/Hour/Minute.

Regarding time message, if you desire to adjust Year/Month/Date and Time setting, please follow the procedures as below:

1. Press "⏸" button at least 5 seconds, Now, you can start year setting when "YEAR" digit blinks.
2. Then you press "▲ MEMO" to forward "YEAR" setting also you can press "PAGE ▼" to backward "YEAR" setting.

3. After finishing Year setting, you can press "TIME" button again.
Now you can adjust Month/Date/Hour and Minute digits blink as you desire.

TROUBLESHOOTING

Phenomenon	Problem	Solution
Nothing is displayed when you push the Power/Start button	Are the batteries exhausted? Are the polarities of batteries in the wrong position?	Replace the batteries with new ones
"Err P" shown on LCD panel	Fail to inflate	1. Please check if the connection on main unit and cuff are loosen. Then connect them. 2. Please check if air tube is broken and change it.
"Err 1" shown on LCD panel	Deflate too fast	Please send it to service station.
"Err 2" shown on LCD panel	1. Movements occur on the arm 2. Interference signal around the BPM	1. Keep still after you push Power/Start button. 2. Please remove interfering sources such as magnetic material or cell phone.
"Err 3" shown on LCD panel	Incorrect measured result	1. You might wear thick jacket and your arm fails to detect pressure values. Please take off your jacket and press BPM again. 2. You might use it incorrectly such as measure pressure value on child

If correct measurement is impossible even after checking the above-mentioned points, consult at the store where you have purchased the unit or the nearest dealer without touching the internal mechanism.

In some very rare cases, there may be error due to the physical condition of the person. In such cases, please consult your doctor.

SPECIFICATIONS

Model: Desktop Arm Blood Pressure Monitor

Measurement: Oscillometric

Measurement Range: Pressure: 20 - 280 mmHg,
Pulse Rate: 40 - 195 beats / min

Accuracy: Pressure: within ± 3 mmHg,
Pulse rate: within $\pm 5\%$ of reading
(subject to the temperature, humidity listed below)

Inflation: Automatic inflation via pump

Pressure Detection: Semi - conductor

Power: 4 X batteries

Operation: Temperature: +10°C to +40°C, Humidity: 30% to 85% RH

Storage: Temperature: -20°C to +60°C, Humidity: 10% to 95% RH

External Dimension: 13 cm (L) x 13 cm (H) x 16 cm (H)

Weight of Main unit: 310 g (without battery), cuff 120 g

Arm Circumference: 22 - 34 cm

HOW TO MAINTAIN AND STORE THE UNIT

1. Use a dry, soft cloth to clean the unit, if desired, or use a cloth lightly dampened with tap water.
2. Do not use alcohol, benzene thinner or other harsh chemicals to clean the device or cuff
3. Before measuring, please wash your hand. Do not wash or wet the arm cuff.
4. Keep away direct sunlight, high humidity or dust for storage.

Note: When the unit will not be used for a longtime, be sure to remove the batteries. Otherwise, the batteries may leak damage the unit.

INFORMATION OF BLOOD PRESSURE

What is blood pressure?

Blood is sent to the arteries with the action of a pump called heart (contraction and dilation). The pressure of blood leaving the heart is called "blood pressure", Blood Pressure pulsate with each beat of the heart. The high blood pressure when the heart contracts are called "Systolic pressure" and the low blood pressure when the heart dilates is called "diastolic pressure", The threshold value for hypertension in adults is defined by the World Health Organization (WHO) as 140/90 mmHg.

Health and blood pressure!

When people reach middle age, the risk of hypertension markedly increase. With aging, the blood vessels age rapidly, Furthermore, because of the obesity and lack of exercise, cholesterol sticks to blood vessels, causing them to lose elasticity, Therefore, watching daily blood pressure help to evaluate our health condition.

Why do we need to monitor blood pressure at home?

By recording the blood pressure values and the measuring conditions such as the measuring time or living state every day, you can know the fluctuation tendency of your blood pressure, which helps control your health. Furthermore, recording of daily blood pressure values is very helpful for your doctor to diagnose.

How to manage your Blood Pressure if it works unsuitable?

If correct measurement is impossible even after checking the above-mentioned points, consult at the store where you have purchased the unit or the nearest dealer without touching the internal mechanism.

In some very rare cases, there may be error due to the physical condition of the person. In such cases, please consult your doctor.

Disposal: The product must not be disposed of along with other domestic waste. The users must dispose of this equipment by bringing it to a specific recycling point for electric and electronic equipment.

For further information on recycling points contact the local authorities, the local recycling center or the shop where the product was purchased. If the equipment is not disposed of correctly, fines or penalties may be applied in accordance with the national legislation and regulations.

GIMA WARRANTY CONDITIONS

Congratulations for purchasing a GIMA product.

This product meets high qualitative standards both as regards the material and the production. The warranty is valid for 12 months from the date of supply of GIMA.

During the period of validity of the warranty, GIMA will repair and/or replace free of charge all the defected parts due to production reasons. Labor costs and personnel traveling expenses and packaging not included.

All components subject to wear are not included in the warranty.


The repair or replacement performed during the warranty period shall not extend the warranty.

The warranty is void in the following cases: repairs performed by unauthorized personnel or with non-original spare parts, defects caused by negligence or incorrect use.

GIMA cannot be held responsible for malfunctioning on electronic devices or software due to outside agents such as: voltage changes, electro-magnetic fields, radio interferences, etc. The warranty is void if the above regulations are not observed and if the serial code (if available) has been removed, cancelled or changed.

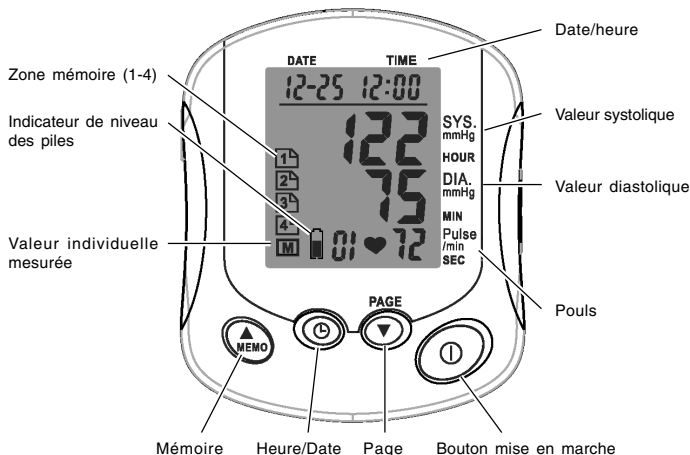The defected products must be returned only to the dealer the product was purchased from. Products sent to GIMA will be rejected.

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