SMARTHEART™

Blood Pressure Monitor with Adult Wrist Cuff

Model # 01-523 series

INSTRUCTION MANUAL

ENGLISH AND SPANISH

Please read this instruction
manual completely before
operating this unit.

STOP!

**PLEASE ENSURE YOU HAVE ALL OF
THE FOLLOWING COMPONENTS
BEFORE USING YOUR DIGITAL
BLOOD PRESSURE MONITOR**


Digital Monitor

Storage Case

Charging Cable

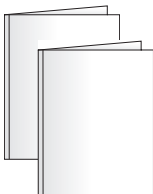Instruction Manual
Quick Start Guide
Blood Pressure Reading Tracking Log
FAQ Information Sheet

**IF YOU ARE MISSING ANY PARTS,
INCLUDING INSERTS OR INSTRUCTION
MANUALS, DO NOT RETURN TO PLACE OF
PURCHASE. CONTACT CUSTOMER CARE AT
866-326-1313.**

INDEX

Care & Safety Information.....	4-7
Introduction	8
Digital Blood Pressure Monitor Features	9
Getting Started	10-11
Charging the Battery	12
Setting Date and Time.....	13-15
Applying Your Blood Pressure Cuff	16-17
Taking Your Blood Pressure Reading	18-21
Interpreting your Results	22-23
Memory Functions	24-25
Care and Maintenance.....	26-27
Device and Label Symbols.....	28
Display Symbols	29
Troubleshooting and Error Codes	30-31
FCC Statement.....	32-33
Electromagnetic Compatibility.....	34-36
Specifications	37
Warranty	38-39
Instrucciones en Español	41-82

**Toll-Free Customer
Care Help Line:
1-866-326-1313**

**Monday – Friday
8:30 a.m. –
4:30 p.m. CST**

SmartHeart™
Manufactured for
Veridian Healthcare
1175 Lakeside Drive
Gurnee, IL 60031

Made in China

#93-1276 12/16
©2016 Veridian Healthcare, LLC

CARE & SAFETY INFORMATION

NOTE: Read all instructions carefully before use. The following basic precautions are needed when using an electrical product.

CAUTION: Failure to read and observe all precautions could result in personal injury or equipment damage.

Improper care or use of your blood pressure monitor may result in injury, damage to the unit or ineffective treatment. Following these instructions will ensure the blood pressure monitor's efficacy and long life.

GENERAL CAUTIONS AND WARNINGS

- Measurement results should be discussed with your physician or healthcare professional; never self-diagnose or attempt treatment as this may be dangerous
- Intended for adult use only; this device is not approved for infant or child use
- Keep out of reach of children
- Do not use the unit if it has any damaged parts, if it has been submersed in water or dropped
- If any abnormality occurs, discontinue use until the unit has been examined and repaired

CARE & SAFETY INFORMATION

- Only use included attachments and parts; do not use attachments from other brands or models as these may not be properly calibrated for use with this device and may result in measurement error
- Do not disassemble the unit or attempt to repair it; substitution of a component different from that supplied may result in measurement error and will void manufacturer warranty
- Always follow local regulations for proper disposal of the monitor, cuff and batteries
- Do not use this device with other diagnostic monitoring equipments
- Take Care to note the blood circulation of the user and discontinue use if prolonged discomfort occurs
- It is inappropriate for people with serious arrhythmia to use this Electronic Sphygmomanometer
- Do not share this device with infected persons where there may be a risk of cross-contamination
- Consult with your physician before use if any of the following conditions are present: arrhythmia, inflammatory diseases, intravascular therapy, arteriovenous shunts or if a mastectomy has been performed

CARE & SAFETY INFORMATION

OPERATING CAUTIONS AND WARNINGS

- Prolonged over-inflation may cause congestion, swelling or bruising in some people
- If you experience any discomfort or an abnormality occurs, stop using the device immediately; to stop the device during operation, press the **POWER** button and the cuff will automatically deflate
- Operation outside of stated operating temperature may result in measurement error or device malfunction; operation environment temperature is: 41°F – 104°F (5°C – 40°C); Humidity: ≤85%RH (non-condensing)
- Never use this unit while operating a vehicle or in the bathtub or shower

STORAGE CAUTIONS AND WARNINGS

- Storage outside of stated storage temperature may result in measurement error or device malfunction; storage environment temperature is: -4°F – 140°F (-20°C – 60°C); Humidity: 10%RH to 93%RH (non-condensing)
- Keep the unit out of reach of small children
- Remove the batteries if the unit will not be used for an extended period of time

CARE & SAFETY INFORMATION

CLEANING CAUTIONS AND WARNINGS

- Never immerse the unit in water to clean as it may damage the unit
- Follow the 'Cleaning and Maintenance' portion of this manual for instruction on how to clean and care for your monitor

INTRODUCTION

It is recommended that you first seek the advice and recommendation of your physician or healthcare professional when using home diagnostic devices, including blood pressure monitors.

The Blood Pressure Monitor is digital monitors intended for use in measuring blood pressure and heartbeat rate with wrist circumference ranging from 5.3" to 8.4" inches (about 13.5 – 21.5 cm). It is intended for adult indoor use only.


Automatic digital blood pressure monitors use the oscillometric method to electronically measure your blood pressure. The monitor detects your blood's movement through the artery in your arm and converts the movements into a digital reading. The oscillometric method does not require a stethoscope, making the monitor ideal for home use.

Blood pressure readings determined with this device are equivalent to measurements obtained by a trained healthcare professional using the cuff/stethoscope auscultation method, within the limits prescribed by the American National Standard for Electronic or Automated Sphygmomanometers.

Thank you for purchasing a SmartHeart Automatic Digital Blood Pressure Monitor. With proper care and use, your monitor will provide you with many years of reliable readings.

**PLEASE READ THIS INSTRUCTION MANUAL
COMPLETELY BEFORE OPERATING THIS UNIT.**

FEATURES

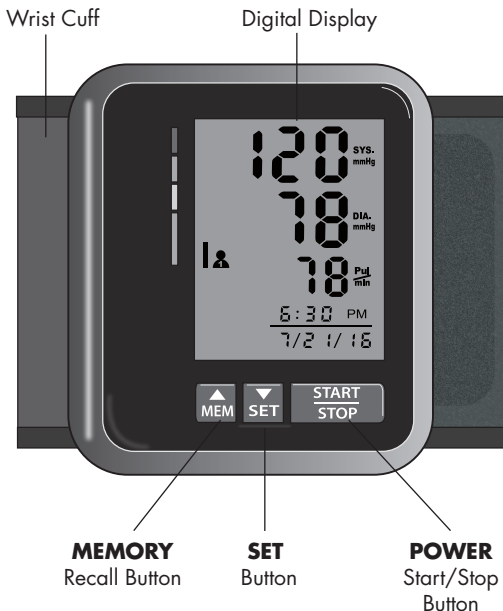DETAILS OF THE DIGITAL DISPLAY SCREEN AND INTERPRETATION OF THE SYMBOLS ARE INCLUDED THROUGHOUT THIS MANUAL.

GETTING STARTED

Before taking a measurement

**NOTE: THIS RECHARGEABLE DEVICE IS
'ASLEEP' FOR SHIPPING AND TRANSPORT.**

To activate your monitor, press and hold the **SET** button for several seconds. The display will illuminate and your monitor is ready for set-up and use.

1. Your readings should only be interpreted by your physician or healthcare professional with access to your individual medical history. Regular use of a home blood pressure monitor will allow you to track and record your readings for discussion with your physician.
2. Conduct your measurement in a quiet place while seated in a relaxed position. Rest for 15 minutes before taking your reading. If an error occurs or you desire to take a second reading, allow 15 minutes between readings for your blood vessels to return to normal.
3. Avoid smoking, eating, taking medication, alcohol, physical activity or any other stressful activity for 30 minutes prior to taking a reading.
4. Always remove any jewelry or constrictive clothing that may interfere with the cuff placement.


GETTING STARTED

5. Keep yourself and the monitor still during measuring; do not talk during the reading.
6. It is recommended that you take your readings at the same time each day to better monitor any indications in your results.
7. Record your daily measurement on the included chart or some other written document to share with your physician.

CHARGING THE BATTERY

This device includes an internal rechargeable battery within the electrical stimulator unit and USB cable. A USB wall charger/adapter is required for charging; not included.

It is necessary to charge the battery when 'Lo' + low battery symbol appears on the display, when the display dims or when the display does not illuminate.

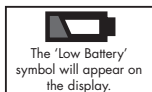
Before first use, the device should be fully charged. Charging time may take up to 90 minutes from depleted battery.

1. Connect the wide end of the USB cable to a wall adapter.
2. Connect the smaller mini USB to the wrist monitor.
3. During charging the digital display will show a battery flashing to indicate the current charge level. When the device is fully charged, the battery symbol will be stop flashing and remain solid on the display.
4. Unplug the charger and stimulator when fully charged; never operate the unit while charging.

SETTING DATE/TIME

The date/time feature is on the bottom right of the display screen. This unit stores the date, time and year allowing for a date/time stamp on all readings in the memory recall function. The **SET**, **MEM** and **START/STOP** buttons are located on the bottom center of the monitor.

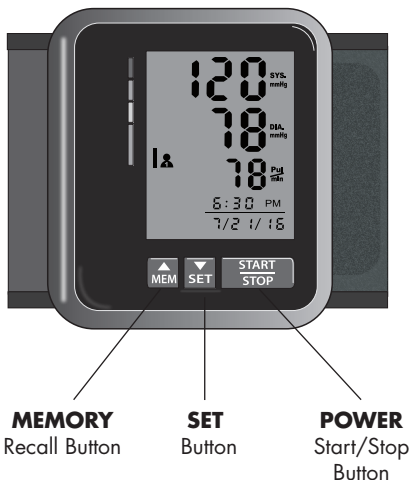SETTING DATE/TIME

HOW TO SET DATE AND TIME

1. When the unit is off (a blank display screen), press and hold the **SET** button for 3-5 seconds to enter clock setting mode.




2. **HOURS/MINUTES**—first the hour will flash on the display. Press the **MEM** button to advance the display to the desired hour, press the **SET** button to confirm your selection and advance to the minute setting. Press the **MEM** button to advance to the correct minute, press **SET** to confirm.

NOTE: the time will show PM when the set time is PM, but the same area will be blank for AM.




3. **MONTH**—next the Month will flash on the display. Press the **MEM** button to advance the display to the desired month, press the **SET** button to confirm your selection and advance to the date settings.
4. **DATE**—repeat the same sequence of step 2 to confirm the MONTH.
5. **YEAR**—repeat the same sequence of step 2 to confirm the YEAR.
6. After confirming YEAR, the LCD will display "dONE" and then shut off.

SETTING DATE/TIME

PLEASE SET DATE AND TIME BEFORE INITIAL USE TO ENSURE YOUR READINGS ARE LABELED AND RECORDED PROPERLY.

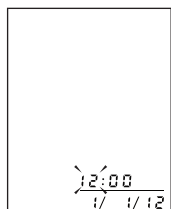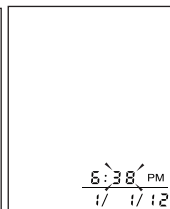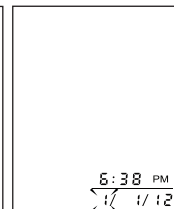		
--	---	---

Hour → Minutes → Month →

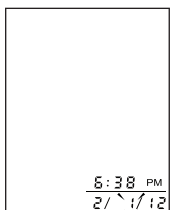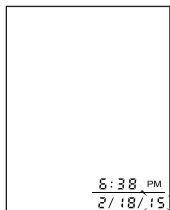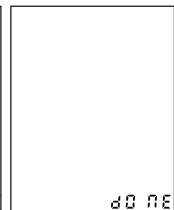		
---	--	--

Date → Year → Complete →

FITTING AND APPLYING YOUR CUFF

PROPER CUFF FIT AND PLACEMENT IS CRITICAL IN OBTAINING ACCURATE BLOOD PRESSURE MEASUREMENTS. PLEASE READ THESE INSTRUCTIONS CAREFULLY AND CONSIDER CONSULTING WITH YOUR PHYSICIAN FOR FURTHER DIRECTION.

Before applying your blood pressure cuff, be sure your wrist is within the appropriate cuff range:


Wrist Circumference Range Suitable for 5.3" – 8.4"
(13.5 cm – 21.5 cm).

Reminder, it is important to avoid smoking, eating, taking medication, alcohol consumption or physical activity 30 minutes prior to taking a reading.

If for any reason you are unable to or cannot not use your left wrist, please modify the instructions for cuff application to your right wrist. Your physician can help you identify which wrist is best for you to take measurements from.

Applying Your Cuff

1. Remove any constrictive clothing or jewelry that may interfere with cuff placement.

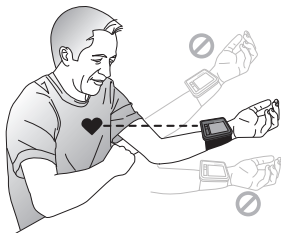FITTING AND APPLYING YOUR CUFF

2. Be seated at a table or desk with your feet flat on the floor.
3. Hold your left arm in front of you with your palm facing upward.
4. Apply the cuff to your left wrist so that the digital display face is positioned on the inside area of your wrist facing you, Fig. 1.
5. Adjust the cuff approximately 0.5" from the edge of the head of the ulna bone, Fig. 2.
6. The cuff should fit comfortably, yet snugly around your wrist.

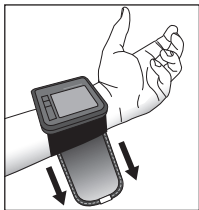
Fig. 1

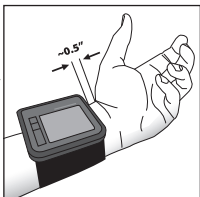
Fig. 2

If you are not comfortable with applying your cuff, please seek the assistance of another member of your household or work with your physician to practice the cuff application. Incorrectly applied cuffs may result in inaccurate readings.

TAKING YOUR BP READING

Please read the preceding portions of this manual prior to taking your first reading.

Reminder: It is important to avoid smoking, eating, taking medication, alcohol consumption or physical activity 30 minutes prior to taking a reading.

Select the desired user, 1 or 2, before measurement.

The user icon displayed represents the memory bank where the measurement result will be stored.

1. When the monitor is OFF, press and hold the **MEM** button. The display will illuminate and the current user, 1 or 2, will flash on the display.
2. Press **MEM** again to toggle between user 1 or 2.
3. Press **SET** to confirm the selection.

Proceed to the measurement steps below.

TAKING YOUR BP READING

Taking the measurement.

1. Rest your elbow on a solid surface with your palm facing upward. Elevate your arm so that the cuff is at the same level as your heart, Fig 1. Be certain that cuff is elevated to the proper level; you may choose to support your hand with a folded towel or book. Relax your left hand. Be seated with your feet flat on the floor and legs uncrossed.

Fig. 1

2. Press and release the **START/STOP** button, the unit will run a self-test, Fig. 2.
3. As the inflation begins, the display will briefly show 0.
4. The cuff will automatically inflate until the reading is detected. You will see the numbers increase on the screen and feel the cuff inflate slowly.

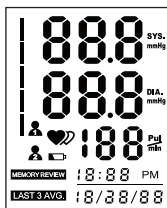
Fig. 2

TAKING YOUR BP READING

YOU CAN STOP THE INFLATION OR DEFLATION PROCESS AT ANY TIME BY PRESSING THE START/STOP BUTTON.

- When the measurement is complete, the cuff will automatically deflate, and your blood pressure measurement and pulse reading results will show on the display screen, Fig 3.
- The Hypertension Indicator will indicate your reading range on the left display in comparison to the color bars on the face of the monitor, Fig 4. See page 23 for information regarding interpreting the Hypertension Indicator results.

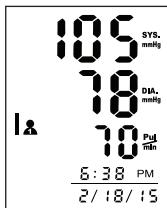
Fig. 3

Fig. 4

TAKING YOUR BP READING

7. If an irregular heartbeat was detected during the reading, the Irregular Heartbeat Detector symbol will appear on the display, Fig 5. 

See page 22 for more information on the Irregular Heartbeat Detector feature.

8. The reading will automatically be stored in memory, up to 60 readings.
9. Press the **START/STOP** button to turn the unit off to conserve energy and battery life. The unit will also automatically shut-off after 1 minute of non-use.



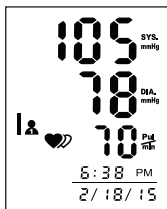
Fig. 5

INTERPRETING YOUR READING

Irregular Heartbeat Detector

An irregular heartbeat is defined as a heartbeat rhythm that is more than 25% slower or 25% faster than the average measurement detected while the monitor is conducting the measurement. When this device detects an irregular rhythm two or more times during the reading, the IHB symbol will appear on the display screen.

This monitor will operate if the icon appears, but the results may be compromised, especially if this symbol appears often. Please consult with your physician or trained healthcare professional for further information regarding an irregular heartbeat and if this symbol appears frequently.

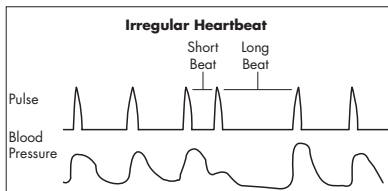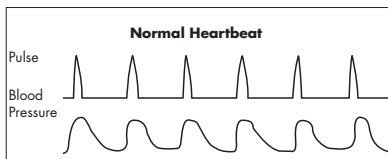INTERPRETING YOUR READING

Hypertension Indicator

This unit features our unique Hypertension Indicator. The color bars on the left side of the monitor display correspond with an icon on the digital display, indicating where the measurement results fall within the American Heart Association standards.

This chart should only be considered as a guideline; always consult with your physician to interpret your individual results. Never adjust your medication(s) dosage, or severely alter your diet or exercise routine without consulting with your physician.

Hypertension Indicator

Category	Systolic (upper #)		Diastolic (lower #)
Hypertensive Crisis	≥180	or	≥110
High Blood Pressure, Stage 2	≥160	or	≥100
High Blood Pressure, Stage 1	140-159	or	90-99
Prehypertension	120-139	or	80-89
Normal	< 120	and	<80

MEMORY FUNCTIONS

This monitor stores and recalls up to 60 readings per user. These measurements can be shared with your physician or trained healthcare professional.

Select the desired user, 1 or 2, for review.

The user icon displayed represents the memory bank where the measurement results are stored.

1. When the monitor is OFF, press and hold the **MEM** button. The display will illuminate and the current user, 1 or 2, will flash on the display.

2. Press **MEM** again to toggle between user 1 or 2.
3. Press **SET** to confirm the selection.

Proceed to the measurement recall steps below.

Recalling Measurements in Memory

1. Press the **MEM** button; the display will show the average of the last three measurements in that memory bank, Fig 1.
2. Press the **MEM** button again to begin scrolling through the saved measurements. Measurements will appear on the display from most current to oldest; the date / time shown will indicate when the measurement was taken.

Fig. 1

MEMORY FUNCTIONS

3. All results for a given measurement will display, including measurement results, pulse rate, Hypertension Indicator and Irregular Heartbeat alert (if applicable), Fig 2.
4. When the number of readings exceeds 60 , the oldest data will be replaced with the new record.
5. Press **POWER** to turn the monitor off at any time during review of the stored measurements.

Fig. 2

Clearing Measurements From Memory

To delete all results in a given bank:

1. When the monitor is off, press and release the **MEM** button to access the memory bank feature.
2. Press and hold both the **MEM** and **SET** buttons for 3 seconds.
3. The display will show 'dEL' and 'dONE', Fig 3.
4. The memory bank is now empty, Fig 4, and the monitor will turn off.

Fig. 3

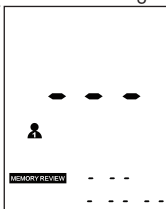
Fig. 4

CARE AND MAINTENANCE

Proper care and maintenance of your blood pressure monitor will help ensure long life for the product; improper use of the monitor may void the included warranty.

Monitor Care and Maintenance

- Do not use any liquids on the monitor; use a soft, dry cloth to clean the monitor as needed
- Storage outside of stated storage temperature may result in measurement error or device malfunction; storage environment temperature is: -4°F – 140°F (-20°C – 60°C); Humidity: 10%RH to 93%RH (non-condensing)
- Avoid high temperatures and direct sunlight
- Remove the batteries if the unit will not be used for an extended period of time
- Do not attempt to disassemble your monitor; disassembling your monitor may void the manufacturers warranty
- Do not subject the monitor to strong shocks; take care not to drop the monitor
- No component can be maintained by user in the monitor. The circuit diagrams, component part lists, descriptions, calibrations instructions, or other information which will assist the user's appropriately qualified technical personnel to repair those parts of equipment which are designated repairable can be supplied

CARE AND MAINTENANCE

Cuff Care and Maintenance

- The cuff may be spot cleaned with a mild-detergent as needed
- Never attempt to iron or dry the cuff; if the cuff becomes wet, allow the cuff to air dry completely before use or storage
- To disinfect the cuff, wipe the inside (skin contact side) of the cuff using a cloth lightly moistened with 70% - 90% ethyl rubbing alcohol and allow to air dry; it is suggested to disinfect the cuff frequently when used in clinical sittings or if the cuff is contaminated in any way.
- It is recommended the performance should be checked every 2 years or after a repair. Please contact customer care at 866-326-1313.
- It is recommended to clean the cuff every 200 uses.
- The monitor can maintain the safety and performance characteristics for a minimum of 10,000 measurements of three years.
- The cuff integrity is maintained after 1,000 open-close cycles of the closure.

DEVICE & LABEL SYMBOLS

These symbols may appear on your device, instructions or packaging and may vary by make and model.

Symbol	Meaning
	Read This Manual—All included manuals should be read prior to first use
	Warning—Symbol indicates a warning, prohibition or mandatory action that mitigates a risk that is not necessarily obvious to the device operator
	Type BF Applied Parts—Indicates that a part of this unit comes in contact with the patient in order to carry out it's intended function; in the case of this device the cuff is the Type BF Applied Part
	Environment Protection—dispose of this product properly; consult with your local recycling ordinances for proper recycling and disposal
	Keep Dry—This device should be kept dry; never submerge the unit or cuff. Consult with the Care and Maintenance section of this manual for information on cleaning your monitor

DESCRIPTION OF DISPLAY SYMBOLS

Symbol	Meaning
SYS	Systolic Blood Pressure—The 'top' blood pressure result
DIA	Diastolic Blood Pressure—The 'bottom' blood pressure result
PUL	Pulse Rate—Pulse rate per minute during the measurement
mmHg	Millimeter of Mercury—Unit of blood pressure measurement for this device
	Low Battery—Battery life indicator
	Pulse Reading—Appears when pulse reading is displayed.
	Irregular Heartbeat Symbol—An irregular heartbeat was detected during the measurement
	User 1/User 2—Indicates which user is active
LAST 3 AVG.	Average—The display shows the average of the last three measurements in memory
MEMORY REVIEW	Memory—The display shows the measurements stored in memory

TROUBLESHOOTING

Problem	Possible Cause	Solution
Display shows an abnormal result	Cuff position was not correct or was not tight enough	Apply the cuff properly and attempt a new measurement
	Body posture was not correct during measurement	Review the 'Taking a Measurement' portion of this manual and attempt a new measurement
	Speaking, excessive movement, angry, excited or nervous anxiety during measurement	Wait a period of time and attempt a new measurement after calm has returned; do not speak or move during measurement
Display is blank when power is on	Battery charge is low	Recharge the device following recharging instructions and attempt a new measurement
Low battery symbol + Lo appear on the display	Battery charge is low	Recharge the device following recharging instructions and attempt a new measurement

ERROR CODES

Code	Meaning	Corrective Action
E1	The cuff is not secure	Re-apply the cuff and attempt new measurement.
E2	The cuff is too tight	Re-apply the cuff and attempt new measurement.
E3	The cuff pressure is in excess	Relax for several minutes and attempt new measurement.
E10 or E11	The monitor detected motion while measuring	Sit comfortable and remain still; attempt a new measurement
E20	The measurement process does not detect the pulse signal.	Re-apply the cuff and attempt new measurement.
E21	The unit measured incorrectly	Relax for a moment and then measure again.
EExx	A calibration error occurred.	Retake the measurement; if the problem persists, contact the Customer Care.

FCC STATEMENT

NOTE:

POTENTIAL FOR RADIO/TELEVISION INTERFERENCE (for U.S.A. only)

This product has been tested and found to comply with the limits for a Class B digital device, pursuant to part 15 of the FCC rules.

These limits are designed to provide reasonable protection against harmful interference in a residential installation. The product generates, uses, and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If the product does cause harmful interference to radio or television reception, which can be determined by turning the product on and off, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna
- Increase the separation between the product and the receiver.
- Connect the product into an outlet on a circuit different from that to which the receiver is connected.

FCC STATEMENT

- Consult the dealer or an experienced radio/TV technician for help.

POTENTIAL FOR RADIO/TELEVISION INTERFERENCE (for Canada only)

This digital apparatus does not exceed the Class B limits for radio noise emissions from digital apparatus as set out in the interference-causing equipment standard entitled "Digital Apparatus", ICES-003 of the Canadian Department of Communications.

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Changes or modifications not expressly approved by the party responsible for compliance could void the user's authority to operate the equipment.

ELECTROMAGNETIC COMPATIBILITY

Table 1 EMISSIONS

For all ME EQUIPMENT and ME SYSTEMS

Guidance and manufacture's declaration - electromagnetic emissions

This monitor is intended for use in the electromagnetic environment specified below. The customer or the user of this monitor should assure that it is used in such an environment.

Phenomenon	Compliance	Electromagnetic Environment Guidance
RF emissions	CISPR 11 Group 1, Class B	Home healthcare environment
Harmonic distortion	IEC 61000-3-2 Class A	Home healthcare environment
Voltage fluctuations and flicker	IEC 61000-3-3 Compliance	Home healthcare environment

ELECTROMAGNETIC COMPATIBILITY

Table 2 — ENCLOSURE PORT

For all ME EQUIPMENT and ME SYSTEMS

Guidance and manufacture's declaration - electromagnetic immunity

This monitor is intended for use in the electromagnetic environment specified below. The customer or the user of this monitor should assure that it is used in such an environment.

Phenomenon	Basic EMC Standard	Immunity test levels Home healthcare environment
Electrostatic Discharge	IEC 61000-4-2	± 8 kV contact ± 2 kV, ± 4 kV, ± 8 kV, ± 15 kV air
Radiated RF EM field	IEC 61000-4-3	10V/m 80MHz-2.7GHz 80% AM at 1kHz
Proximity fields from RF wireless communications equipment	IEC 61000-4-3	Refer to table 3
Rated power frequency magnetic fields	IEC 61000-4-8	30A/m 50Hz or 60Hz

ELECTROMAGNETIC COMPATIBILITY

Table 3 — PROXIMITY FIELDS FROM RF WIRELESS COMMUNICATIONS EQUIPMENT

For ME EQUIPMENT and ME SYSTEMS that are not LIFE-SUPPORTING

Guidance and manufacturer's declaration - electromagnetic immunity

This monitor is intended for use in the electromagnetic environment specified below. The customer or the user of this monitor should assure that it is used in such an environment.

Test frequency (MHz)	Band (MHz)	Immunity test levels Professional healthcare facility environment
385	380-390	Pulse modulation 18Hz, 27V/m
450	430-470	FM, ± 5 kHz deviation, 1kHz sine, 28V/m
710	704-787	Pulse modulation 217Hz, 9V/m
745		
780		
810		
870	800-960	Pulse modulation 18Hz, 28V/m
930		
1720		
1845	1700-1990	Pulse modulation 217Hz, 28V/m
1970		
2450	2400-2570	Pulse modulation 217Hz, 28V/m
5240	5100-5800	Pulse modulation 217Hz, 9V/m
5500		
5785		

PRODUCT SPECIFICATIONS

Name	Blood Pressure Monitor with Adult Wrist Cuff
Model Number	01-523 series
Display System	Digital display/LCD
Measuring Method	Oscillometric Method, Automatic Air Inflation and Measurement
Power Source	3.7V 420mAH Built-in rechargeable lithium-ion battery, 5V  1A USB cable
Measuring Range	Cuff pressure: 0-300 mmHg; Measurement Pressure: 40-230 mmHg; Pulse: 40-199 beats/minute
Accuracy	Within operating environment 41°F – 104°F (5°C – 40°C): Pressure: ±3 mmHg; Pulse: ±5% of readings
Pressurization	Automatic pressurization by pump
Deflation	Automatic pressure release valve
Memory	Built-in memory with 2 banks, 60 measurements per bank
Automatic Power-Off	Approximately 1 minute after last button operation
Battery Life	Approximately 200 measurements
Operation Environment	Temperature 41°F – 104°F (5°C – 40°C); Humidity ≤85% RH (non-condensing)
Storage Environment	Temperature -4°F – 140°F (-20°C – 60°C); Humidity 10%RH to 93%RH
Monitor Dimensions	3.14" x 2.85" x 0.5"
Monitor Weight	3.8 oz. (without batteries)
Wrist Circumference	Suitable for 5.3" – 8.4" (13.5 – 21.5 cm)
Accessories	USB cable, instruction manual, storage case, Quick Start Guide, Blood Pressure Log

Specifications are subject to change without notice

LIFETIME LIMITED WARRANTY

Congratulations on your purchase of a Digital Blood Pressure Monitor. Your Digital Blood Pressure Monitor is covered by the following limited warranty commencing upon the date of purchase, and subject to the following terms and conditions:

The Warranter warrants that its Digital Blood Pressure Monitor will be free from defects in materials and workmanship under normal consumer usage for the life of the product for the original purchaser.

Periodic maintenance, repair and replacement of parts due to normal wear and tear are excluded from coverage. Defects or damage that result from: (a) improper operation, storage, misuse or abuse, accident or neglect, such as physical damage (cracks, scratches, etc.) to the surface of the product resulting from misuse; (b) contact with liquid, water, rain, extreme humidity or heavy perspiration, sand, dirt or the like, extreme heat, or food; (c) use of the Digital Blood Pressure Monitor for commercial purposes or subjecting the Digital Blood Pressure Monitor to abnormal usage or conditions; or (d) other acts which are not the fault of the Warranter, are excluded from coverage. This warranty does not cover batteries or other power sources that may be provided with, or used with the Digital Blood Pressure Monitor.

If the Digital Blood Pressure Monitor fails to conform to this limited warranty, return the Digital Blood Pressure Monitor postage prepaid to: **Attn: Repair Department, 1175 Lakeside Drive, Gurnee, IL, 60031**. When returning a product, please also include: (i) a copy of your receipt, bill of sale or other comparable proof of purchase; (ii) a written description of the problem; and (iii) your name, address and telephone number. Carefully package the product to avoid any damage that may occur while in transit; shipping insurance with returned receipt is recommended. At our option, the Warranter will repair or replace the unit found to be defective in materials or workmanship under normal consumer usage. The purchaser will be notified of any additional repairs required prior to completing the repair, and will

LIFETIME LIMITED WARRANTY

be responsible for parts charges, if any, and repair charges not covered by this limited warranty.

EXCEPT AS PROVIDED FOR IN THIS LIMITED WARRANTY, ALL EXPRESS AND IMPLIED WARRANTIES AND CONDITIONS ARE DISCLAIMED, INCLUDING WITHOUT LIMITATION THE IMPLIED WARRANTIES OF MERCHANTABILITY AND FITNESS FOR A PARTICULAR PURPOSE. THE REPAIR OR REPLACEMENT AS PROVIDED UNDER THIS LIMITED WARRANTY IS THE EXCLUSIVE REMEDY OF THE CONSUMER, AND IS PROVIDED IN LIEU OF ALL OTHER WARRANTIES, EXPRESS OR IMPLIED. IN NO EVENT SHALL THE WARRANTER BE LIABLE, WHETHER IN CONTRACT OR TORT (INCLUDING NEGLIGENCE) FOR DAMAGES IN EXCESS OF THE PURCHASE PRICE OF THE PRODUCT, OR FOR ANY INDIRECT, INCIDENTAL, SPECIAL OR CONSEQUENTIAL DAMAGES OF ANY KIND, OR FOR DAMAGES TO, OR LOSS OF, OTHER PROPERTY OR EQUIPMENT OR PERSONAL INJURIES TO THE FULL EXTENT THESE DAMAGES MAY BE DISCLAIMED BY LAW.

Some states and jurisdictions do not allow the limitation or exclusion of incidental or consequential damages, or limitation on the length of an implied warranty, so the above limitations or exclusions may not apply to you. This warranty gives you specific legal rights, and you may also have other rights that vary from state to state or from one jurisdiction to another.

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