
# ***Sterling Series Floor-Type Clock Aneroid Sphygmomanometer***

Model 02-131

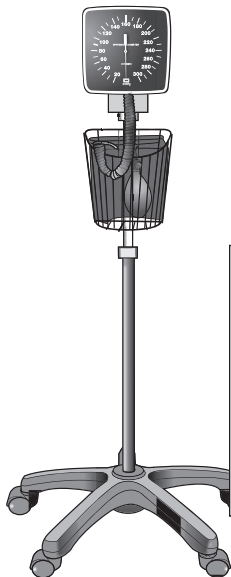Optional wall-mount  
hardware included

**INSTRUCTION MANUAL • ENGLISH & ESPAÑOL**

---

# INDEX

---

|                                         |       |
|-----------------------------------------|-------|
| Introduction .....                      | 3     |
| Parts and Mounting Hardware .....       | 4-5   |
| Floor Stand Assembly Instructions ..... | 6-7   |
| Wall Mount Instructions .....           | 8     |
| Clock and Cuff Assembly .....           | 9     |
| Gauge Adjustment Feature .....          | 10    |
| Arm Cuff Size and Selection .....       | 11    |
| Care and Maintenance .....              | 11    |
| Instrucciones en Español.....           | 13-24 |

---

**Toll-Free Customer Care Help Line: 1-866-326-1313**  
**Monday – Friday 8:30 a.m. – 4:30 p.m. CST**

Distributed for:  
Veridian Healthcare  
1465 S. Lakeside Drive  
Waukegan, IL 60085  
[www.veridianhealthcare.com](http://www.veridianhealthcare.com)



Made in China  
#93-1076 02/10  
©2010 Veridian Healthcare, LLC

# INTRODUCTION

---


Thank you for purchasing a Veridian Healthcare Sterling Series Clock Aneroid Sphygmomanometer. With proper care and use, your sphygmomanometer will provide you with many years of reliable readings.

This clock aneroid sphygmomanometer includes hardware for both a mobile floor stand set-up and wall mount installation. The following instructions include detailed illustrations for both floor stand and wall-clock assembly.

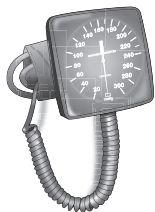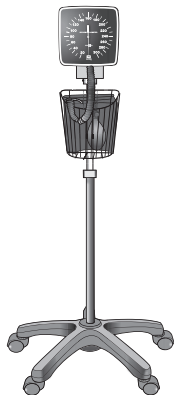**PLEASE READ THIS INSTRUCTION MANUAL  
COMPLETELY BEFORE ASSEMBLING AND  
OPERATING THIS UNIT.**


## INCLUDED PARTS AND HARDWARE
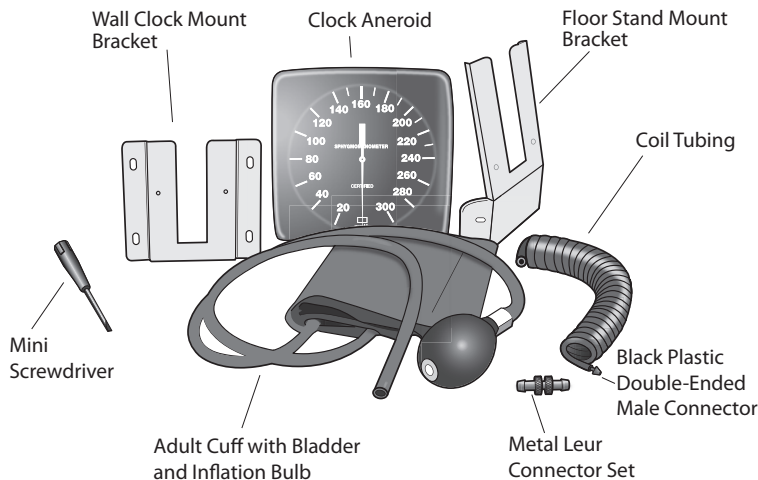
### **STOP!**

**PLEASE ENSURE YOU HAVE ALL OF THE LISTED COMPONENTS BEFORE ATTEMPTING TO ASSEMBLE YOUR ITEM. IF YOU ARE MISSING ANY PARTS, DO NOT RETURN TO PLACE OF PURCHASE. CONTACT CUSTOMER CARE AT 866-326-1313.**

# INCLUDED PARTS AND HARDWARE
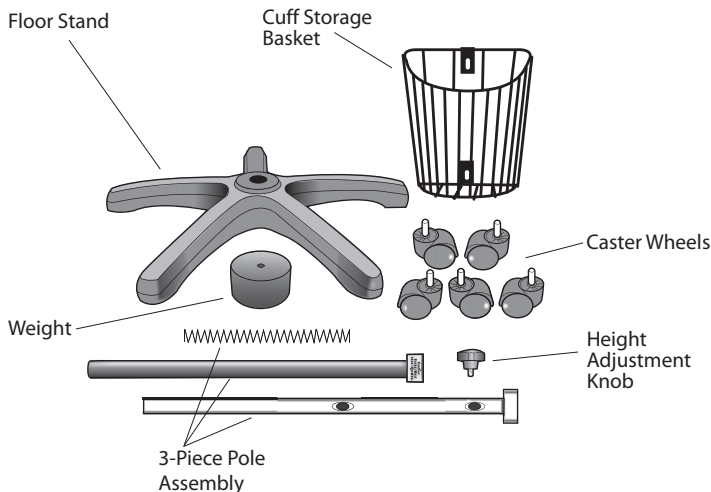
## COMPONENTS:

Wall Clock  
Floor Stand Mount Bracket  
Wall Clock Mount Bracket  
Adult Size Nylon Cuff with Bladder and Inflation Bulb  
Black Plastic Double-Ended Male Connector (shown attached to coil tubing)  
Metal Leur Connector Set  
Coil Tubing  
Mini Screwdriver  
Floor Stand  
Five (5) Caster Wheels  
Metal Basket  
3-piece Pole Assembly (black pole, silver pole, spring)  
Height Adjustment Knob  
7.5 lb. Weight

## HARDWARE:

### Floor Stand Hardware

One (1) 4.8/M8 Screw Bolt with 2 Washers  
Two (2) 3/8" Screws with Integrated Washers + Two (2) Washers  
Two (2) 1/2" Screws

### Wall Mount Hardware

Four (4) 1" Screws  
Four (4) Drywall Anchors  
Two (2) 1/4" Lock Screws

## TOOLS REQUIRED:

Phillips-Head Screwdriver  
Socket Wrench with 9/16" Attachment  
Power Drill (for wall installation option only)

# FLOOR STAND ASSEMBLY INSTRUCTIONS

## ASSEMBLY OF THE STAND/BASE REQUIRES A PHILLIPS-HEAD SCREWDRIVER AND SOCKET WRENCH WITH 9/16" OR 14MM ATTACHMENT.

1. Insert the 5 caster wheels into the black plastic base. The casters will 'click' into place, Fig 1. Gently pull on each one to ensure they are completely attached. Turn the base over onto the wheels.
2. Insert the black pole into the base; smaller end down, open end up, Fig 2.
3. Turn the pole and base over to attach the weight to the center of the underside of the base. Thread the two washers onto the bolt, and insert the bolt through the weight and into the stand base, securing with the socket wrench, Fig 3.
4. Turn the stand over so that the fixture is resting upright on the caster wheels. Insert the spring into the black bottom tube, Fig 4.
5. Insert the silver top portion of the pole into the black bottom portion, Fig 5.
6. Attach the height adjustment knob to stand, turning clockwise to tighten, counterclockwise to loosen, Fig 6.

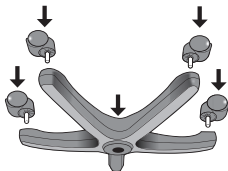
Fig. 1

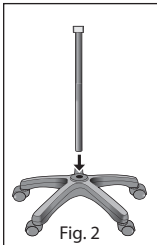
Fig. 2

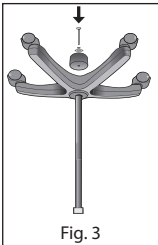
Fig. 3

Fig. 4

Fig. 5

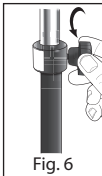
Fig. 6

# FLOOR STAND ASSEMBLY INSTRUCTIONS

7. Using the two (2) 3/8" screws with washers, attach the black wire basket to the upper pole with the phillips-head screwdriver, Fig 7.

**NOTE:** To prevent tipping, align the front of the basket with the caution warning label on the lower pole

8. Attach the floor stand mount bracket to the top of the stand assembly using the two (2) 1/2" screws and the phillips-head screwdriver, Fig 8.
9. Carefully slide the clock fixture onto the bracket, Fig 9. Take care to properly align the brackets as the clock is not further secured with screws or other hardware.
10. Adjust the stand to the desired height using the height adjustment knob. The stand can be adjusted from 45"–55". Never adjust the fixture higher than 55" to prevent tipping or damage to the unit.

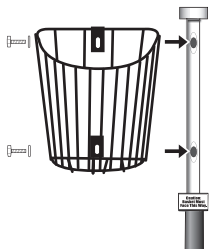
Fig. 7

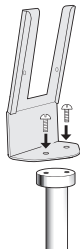
Fig. 8

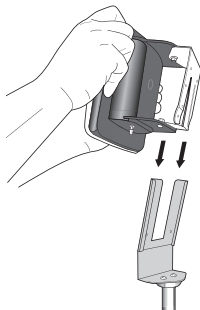
Fig. 9

# WALL MOUNT ASSEMBLY INSTRUCTIONS

**ASSEMBLY OF THE WALL MOUNT OPTION REQUIRES A PHILLIPS-HEAD SCREWDRIVER AND POWER DRILL.**

**KEEP IN MIND THE DESIRED LOCATION OF PATIENT CARE WHEN SELECTING THE LOCATION TO MOUNT THE CLOCK ANEROID FIXTURE.**

1. Using the wall mount bracket as template, use a pencil to mark the screw hole location directly onto the wall.
2. Drill four pilot holes into the wall using a drill bit smaller than the width of the drywall anchor over each pencil mark.
3. Insert the drywall anchors into the four holes until they are flush with the wall; it may be necessary to use a hammer or mallet to lightly tap the anchors into place.
4. Secure the wall mount bracket to the drywall anchors using the four 1" Phillips-Head screws; do not over tighten, Fig 1.
5. Slide the wall clock onto the wall mount bracket, Fig 2.
6. Secure the wall clock onto the bracket using the two 1/4" lock screws, Fig 3.

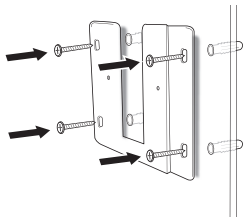
Fig. 1

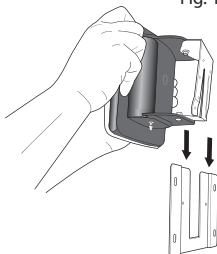
Fig. 2


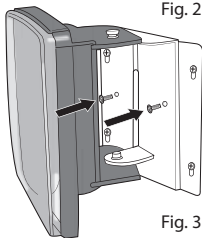
Fig. 3


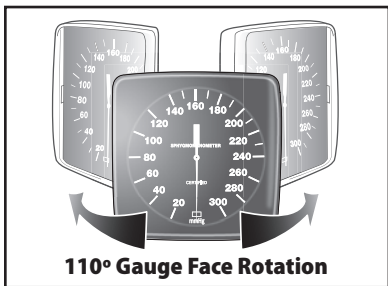
---

## CUFF, BLADDER AND CLOCK ASSEMBLY


---

After the clock is secured to the either the floor stand or wall mount assembly, attach the coil tubing, cuff and bulb as illustrated.

Two connector options are included for user preference: Black Plastic Double-Ended Male Connector and Metal Leur Connector Set (shown in illustration)

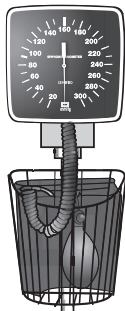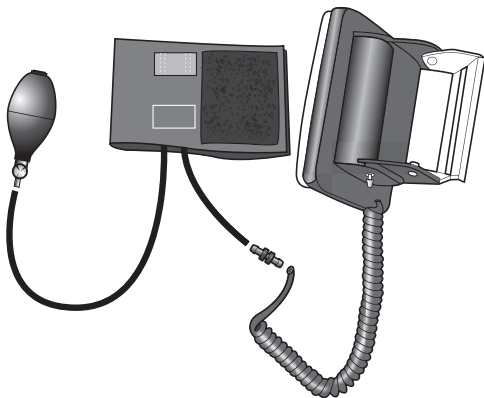Floor Stand – Store the arm cuff in the attached wire basket when not in use.

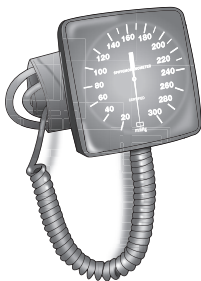Wall Mount – Store the arm cuff in the space between the Wall Clock and the Wall Mount Bracket when not in use.

# GAUGE ADJUSTMENT FEATURE

This unit includes a small screwdriver, which allows the gauge to be adjusted and set to zero when the needle is not within the zero square position.

*Gauge out of zero square position*

*Gauge properly set to zero square position*

*The mini screwdriver is specifically sized for this task. Take care to store the screwdriver in an easily located place.*

1. Remove the cuff and tubing from the air stem on the underside of the gauge, Fig 1.
2. Carefully insert the small screwdriver into the stem until resistance is felt; slowly turn the screwdriver, you will feel the screwdriver connect with the grooves in the adjustment mechanism in the clock, Fig 2.
3. Slowly turn the screwdriver again while observing the needle on the gauge face. Turn until the needle has returned to the center of the zero square.
4. Remove the screwdriver and reattach the coil tubing to the air stem.

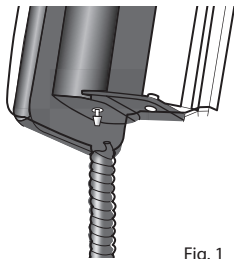
Fig. 1

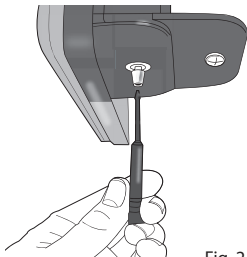
Fig. 2

Resetting the needle to the zero position does not guarantee that the gauge is properly calibrated. Calibration can only be confirmed when a gauge is calibrated against a calibrated mercury column, however, wall-clock gauges are commonly accurately calibrated when the needle is at zero.

---

## ARM CUFF SIZE AND SELECTION

---

This product includes an Adult sized nylon cuff.  
Additional cuff sizes are available through your dealer.

| Size        | Dimensions<br>(Inches) | Arm Girth<br>(Inches) | Dimensions<br>(Centimeters) | Arm Girth<br>(Centimeters) |
|-------------|------------------------|-----------------------|-----------------------------|----------------------------|
| Infant      | 3.25"W x 10"L          | 5.5" – 7.75"          | 8.25cm W x 25.4cm L         | 14cm – 19.7cm              |
| Child       | 4.25"W x 13.5"L        | 7.75" – 11.375"       | 10.8cm W x 34.3cm L         | 19.7cm – 28.9cm            |
| Adult       | 5.5"W x 21"L           | 11" – 16.375"         | 14cm W x 53.3cm L           | 27.9cm – 41.6cm            |
| Large Adult | 7"W x 25.5"L           | 13" – 20"             | 17.8cm W x 64.8cm L         | 33cm – 50.8cm              |
| Thigh       | 8.25"W x 30.5"L        | 16.125" – 24.125"     | 20.9cm W x 77.5cm L         | 41cm – 61.2cm              |

---

## CARE AND MAINTENANCE

---

**Clock Gauge** – To clean, wipe clock face, housing and fixture with a slightly damp cloth.



**Blood Pressure Cuff** – Always remove the bladder prior to cleaning; cuffs may be washed with mild detergent and water or wiped with a damp cloth.

Never use steam or liquid disinfectant to clean any parts of this unit. Do not autoclave any parts of this unit.

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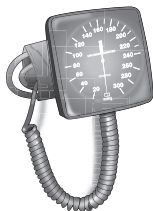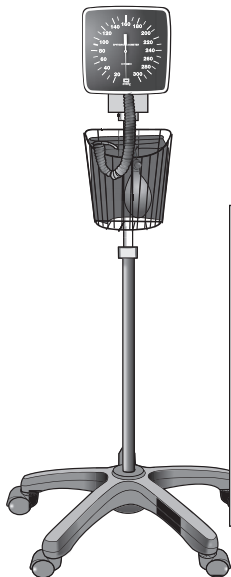Distributed for:  
Veridian Healthcare  
1465 S. Lakeside Drive  
Waukegan, IL 60085  
[www.veridianhealthcare.com](http://www.veridianhealthcare.com)  
Made in China  
#93-1076 02/10  
©2010 Veridian Healthcare, LLC

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