Norditalia

IT MANUALE D'USO

BP-1000

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TABLE OF CONTENTS

Introduction 4
Symbols Used 4
Important Safety Notes 5
Appliance Description 6
Irregular Heartbeat (IHB) Indicator 8
 Irregular Heartbeat Symbol ((♥)) 8
 Arrhythmia 8
Before Taking A Measurement 9
Battery Installation / Replacement 10
Settings 10
 1. Date and time 10
 2. Save Settings 11
Proper Way of Measurement 11
Take a Measurement 13
Classification of Blood Pressure 13
Use the Memory Function 14
 To read the memory values 14
 To delete the memory value 14
About Blood Pressure 15
 Blood Circulation 15
 Health and Blood Pressure 15
 Hypertension 15
 Hypotension 16
Exceptional Situations 16
 Error Indicators 16
 Fault Conditions 16
Care 17
Maintenance 17
Specifications 18
Warranty Conditions and Customer Care 19

0900378003 ed. 09/2015

Introduction

Thank you for purchasing Norditalia Fully Automatic Arm Digital Blood Pressure Monitor.

The appliance is based onto the oscillometric principle for blood pressure measurement.

Please read this instruction manual thoroughly before using the appliance. For specific information about your own blood pressure, consult your doctor.

Symbols Used


SYMBOLS USED IN THIS INSTRUCTION MANUAL AND ONTO THE APPLIANCE	
	WARNING: read carefully the instruction manual
	Type BF applied part
	Manufacturing Date
CE 0051	Compliant to european directive 93/42/EEC.
	Manufacturer for 93/42/EEC directive
LOT	Lot number
	Dispose of the appliance and its accessories in compliance with the regulations in force. Do not dispose of this appliance in the environment!
	Keep dry

⚠ Important Safety Notes

- Read all the instructions carefully before using the appliance.
- Do not use the appliance for any purpose other than measuring blood pressure on adult. Do not use onto neonatal baby.
- Always consult your doctor for specific information about your blood pressure. Self-diagnosis and treatment using measured results may be dangerous.
- Never change your therapy without consulting your doctor.
- Do not measure your blood pressure after bathing, drinking alcohol, smoking,eating or physical exercise.
- Consult your doctor during pregnancy, arrhythmia and arteriosclerosis. People with severe blood flow problems, or blood disorders, should consult a doctor before using the appliance. Cuff inflation can cause internal bleeding.
- Do not leave the appliance unattended with infants or persons who cannot express their consent.
- If battery fluid should get on your skin or eyes, immediately rinse with plenty of clean water. Consult a doctor immediately.
- Use only alkaline batteries with this appliance. Respect the polarity. Remove immediately depleted batteries. Do not use old and new batteries together. Remove batteries if the appliance will not be used for three months or more.
- Do not disassemble or attempt to repair the appliance or components.
- Do not shock or drop the appliance.
- Do not wet the appliance and accessories including cuff.
- Do not inflate the arm cuff when it is not wrapped around your arm.
- Do not use a mobile phone, or other devices that emit electromagnetic fields, near the appliance. This may result in incorrect operation of the appliance.
- Do not use the appliance where flammable gas (such as anesthetic gas, oxygen or hydrogen) or flammable liquid (such as alcohol) are present
- Use only the original AC adapter. When using AC adapter (optional) never plug in or unplug the power cord from the electric outlet with wet hands. Plug the AC adapter into the appropriate voltage outlet. Do not use the AC adapter if the unit or the power cord is damaged.

Appliance Description


Body & Cuff

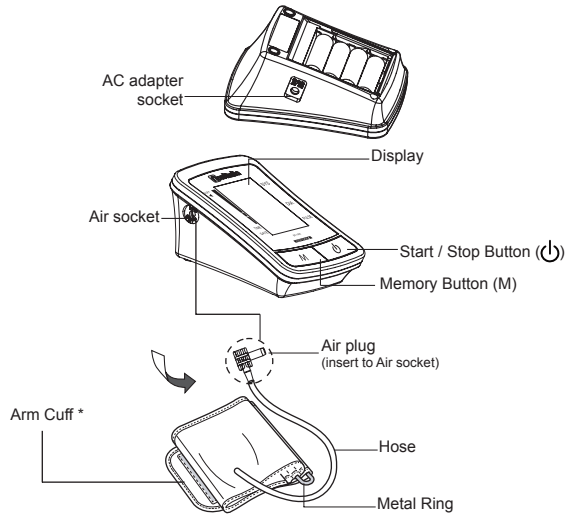*** BP-1000**: Standard Arm Cuff: fit for 22-36 cm range of arm circumference.

To guarantee the performances, functionalities and safety of the appliance, only use Norditalia original accessories and spare parts.

BP-1000 Display

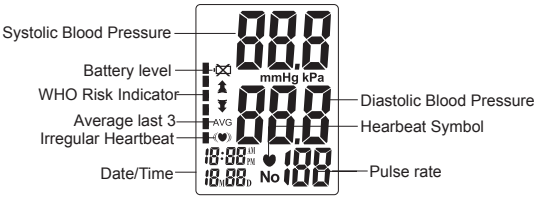
Accessories


4 x AA (LR6) Alkaline Batteries


Storage Bag

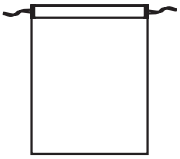AC adapter (optional) – available on request

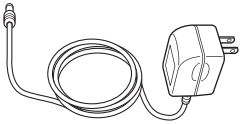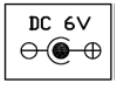Irregular Heartbeat (IHB) Indicator

Irregular Heartbeat Symbol ((♥))

The appliance can detect irregular heartbeats.

When the appliance detects an irregular rhythm twice or more times during the measurement, the irregular heartbeat Symbol will appear on the display with the measurement values.

An irregular heartbeat rhythm is defined as a rhythm that varies by more than 25% from the average rhythm detected during the measurement.

If the Irregular Heartbeat Symbol displays often with your measurement results, we recommend you consult your doctor.

Arrhythmia

Arrhythmia is a condition where the heartbeat rhythm is abnormal because of defects in the electrical stimulation of the heart.

Typical symptoms are: skipped heartbeats, premature contraction, too rapid pulse (tachycardia - atrial fibrillation) or too slow pulse (bradycardia).

A frequent appearance of the irregular heartbeat symbol could indicate arrhythmia but this can only be determined by your doctor through a specific test.


Self-diagnosis and self treatment are always dangerous, consult your doctor.

Before Taking A Measurement


1. Avoid eating,drinking alcohol or coffee, smoking, exercising and bathing for at least 30 minutes before taking a measurement.
2. Stress raises blood pressure. Avoid taking measurements during stressful times.
3. The cuff can be applied to your left or right arm. Always measure your blood pressure on the same arm. If there was high differences between the two arms ask you doctor which arm you should use.
4. Measurements should be taken in a quiet place.
5. Position the appliance at heart level throughout the measurement.
6. Remain still and do not talk during the measurement.
7. A single measurement does not provide an accurate indication of your true blood pressure. You need to take several readings over a period of time. Usually it is suggested to take 3 measurements and calculate the average. Please note the appliance can calculate automatically the average of last 3 measurements (AVG). Please see “Use the memory function”.
8. Your blood pressure varies during the day. Try to measure your blood pressure at the same time each day for consistency (for example within one hour after waking up).
9. It is normal for two measurements close in time to give even very different values.
10. Wait 2 or 3 minutes between measurements. The wait time allows the arteries to return to normal conditions.
11. Measurement taken by your doctor can be different from the measurements you get at home as the situation is different.

Battery Installation / Replacement


Remove the battery cover

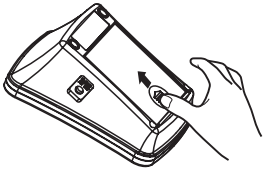Install four “AA / LR6” size batteries so the + (positive) and – (negative) polarities match the polarities indicated in the battery compartment.

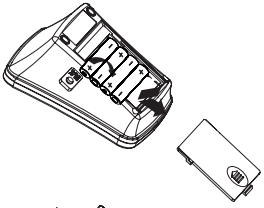Replace the battery cover.
NOTE: stored measurements will not be deleted during battery replacement

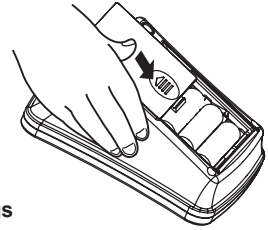Settings

1. Date and time

With power off, press “⏻” button about 3 seconds to actuate system setting.

Set the month first by adjusting the “M “ button. Press ”⏻ ” button again to confirm current month. Continue setting the day, hour and minute in the same fashion.


2. Save Settings

While in any setting mode, press "⏻" button about 3 seconds to turn the unit off. All information will be saved.


The unit will automatically turn off after 3 minutes if not used.

Proper Way of Measurement


1. Make sure the Air plug is inserted properly in the Air socket
2. Remove all clothes from your arm allowing the cuff to fit directly on the skin.

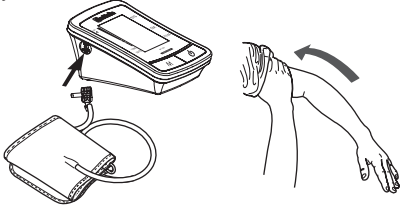
3. Insert the final part of the cuff in the metal ring and keep the hose outward.
4. Put your arm in the cuff, then pull it up till the upper part of your arm.

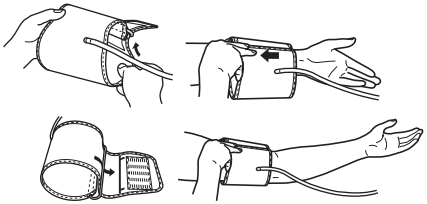5. The hose should be in line with your inner arm. The bottom border of the cuff should be approximately 1-2 cm above your elbow cavity.

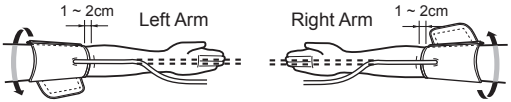
6. Wrap the cuff tightly around your arm by using the Velcro strip.
7. Sit with your feet flat on the floor. Place your arm on a table so the arm cuff will be at heart level during the measurement.
8. Relax your arm, keep the palm upward and the fingers curving naturally.

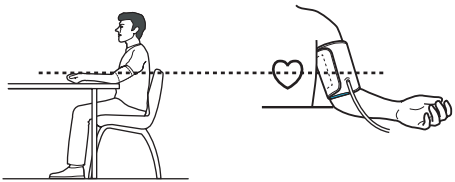Take a Measurement

1. Press the "⏻" button. All display symbols appear.
2. The cuff will inflate to a suitable pressure (measurement starts during inflation to prevent excessive inflating).
3. The cuff will then deflate slowly and pressure and pulse rate will be measured. When an heartbeat is detected the heartbeat symbol will flash.

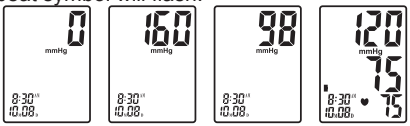1-Start 2-Inflating 3-Deflating 4-Result

4. The blood pressure and pulse rate are displayed when the measurement is completed.
5. The cuff is deflated automatically, and all of the measurement results are stored in the memory with date and time. The IHB symbol ((♥)) will be displayed if irregular heartbeat is detected.
6. Press the "⏻" button to turn off the unit.


The appliance will automatically turn off after 3 minutes if not used.

NOTE: To stop the inflation or measurement, press the "⏻" button anytime. The unit will stop inflating and then release and turn off.

Classification of Blood Pressure

After each measurement is completed, the LCD display will show your position automatically on the six segments of the bar indicator which corresponds to World Health Organization (WHO) Risk Indicator.

Reference Material: Journal of Hypertension 1999, Vol 17 No.2

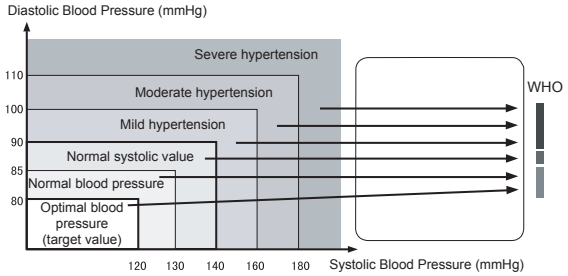Note:

When a person's systolic and diastolic pressures fall into different categories, the higher category should apply.


Use the Memory Function

After every measurement the appliance memorizes the measured blood pressure and pulse rate with date and time. It can store up to 120 sets of measurements

To read the memory values

With power off, press the "M" button to activate screen display. The device will display the average of the last 3 readings. The "AVG" symbol will appear.

Press "M" again to cycle through memories.

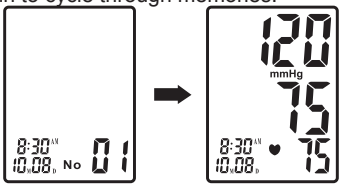To delete the memory value

While in memory check mode. Press and hold the "⏻" button for approximately 3 seconds to delete all memory records. The monitor will beep indicating successful deletion and then transfer into testing mode.


About Blood Pressure

Blood Circulation

The blood circulation is responsible for supplying the body with oxygen. Blood pressure is the pressure exerted on blood vessels walls. The systolic blood pressure value (top value) represents the blood pressure produced by contraction of the heart muscle. The diastolic blood pressure value (lower value) represents the blood pressure produced by relaxation of the heart muscle.

Health and Blood Pressure

The incidence of hypertension increases with age. In addition, a lack of exercise, excess of body fat and high levels of cholesterol (LDL), which sticks to the inside in blood vessels, reduces elasticity of these vessels. Hypertension accelerates arterial sclerosis which can lead to vary serious conditions such as stroke and myocardial infaction. For these reasons it is very important to know wheter the blood pressure is within a healty range. Blood pressure fluctuates from minute to minute, throughout the day. Therefore it is essential to take regular measurements to help you identify an average blood pressure.

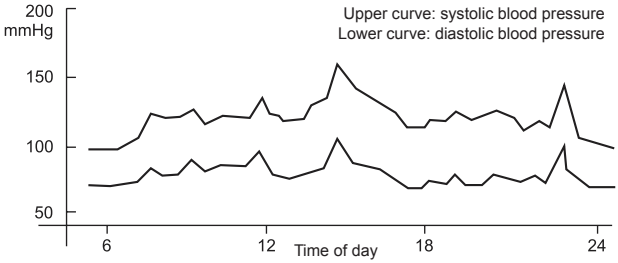Hypertension

High blood pressure can go unnoticed for a long time, since it does not produce noticeable symptoms. If your blood pressure falls in an “hypertension” area in repeated measurements over

several days, you should consult your doctor for detailed medical examination. Self- iagnosis and treatment using measured results may be dangerous. Never initiate, change or suspend a therapy without consulting your doctor.

Hypotension

The definition of hypotension is not universally accepted. Nevertheless, a systolic pressure lower than 100 mmHg is assumed as hypotensive.

Exceptional Situations

Error Indicators

The following symbols will appear on the display in case of misfunction.

Problem	Possible cause	Remedy
	Cuff fails to inflate properly	Make sure hose is properly fastened to cuff and monitor unit
	Improper operation	Read user manual carefully and re-test properly.
	Pressurization is over 300mmHg	Read user manual carefully and re-test properly.
	Low battery power	Replace all of the worn batteries with new ones.

Fault Conditions

Problem	Causes and Solutions
No power. No display appears	Check the battery installation for proper placement of the batteries polarities. If necessary replace batteries.
Measurement values appear too high or too low.	Blood pressure varies constantly. Many factors including stress, time of day, how you wrap the cuff, may affect your blood pressure. Review the sections “Proper Way of Measurement”. Repeat the measurement

Care

Keep the unit in the storage case when not in use.
Clean the unit with a soft dry cloth. Do not use any abrasive or volatile cleaners.
Never immerse the unit or any components in water.


Avoid subjecting the monitor to strong shocks or vibrations, such as dropping the unit on the floor.
Use the unit consistent with the instructions provided in this manual. Use only authorized parts and accessories.

Maintenance

Do not clean the body and cuff with naphtha, thinner or gasoline, etc.
Do not wet the cuff or attempt to clean the cuff with water.
Store the unit in a clean and dry location. Do not subject the unit to extreme hot or cold temperature, humidity and direct sunlight.
Remove the batteries if the unit will not be used for 3 months or longer. Always replace all of the batteries with new ones at the same time.

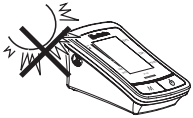Specifications

Model	BP-1000	
Display	LCD Digital Display	
Measuring principle	Oscillometric method	
Measurement position	Upper arm	
Measurement Range	**Pressure**	0-300 mmHg
	Pulse	30-180 Pulses/min
Accuracy	**Pressure**	±3 mmHg
	Pulse	±5% of reading
Inflation and deflation	Automatic	
Memory function	120 memories	
Power source	4 x AA/LR6 Alkaline Batteries	
Automatic turn off	In 3 minutes	
Main unit weight	Approx. 265g. (batteries not included)	
Wrist Cuff	Standard Adult Cuff, suitable for 22-36 cm arm circumference.	
Operating Environment	**Temperature**	5 ÷ 40°C
	Relative Humidity	15 ÷ 85%RH Noncondensing
	Atmospheric Pressure	700 ÷ 1060 hPa

Specifications are subject to change without notice

Warranty Conditions and Customer Care

The product is guaranteed for 2 years from the purchase unless differently required by local regulations against manufacturer's defect under normal, household use with respect of present manual instructions. The warranty is valid only if the device is sent back to the supplier together with the receipt proving the date of purchase and the warranty card duly filled. The warranty does not cover consumables, expiring components and accessories other than the device itself (ex: batteries, cuff, etc...).

For assistance in obtaining accessories or about the use of the device, please fill the form available in the "assistance" section on the website www.norditalia.biz or contact us using the referrals indicated in the manual. If necessary contact your dealer. Periodically (every 2 years), or in case of doubt about the accuracy of the appliance, compare your results with a different appliance. In case of systematic differences higher than 30 mmHg contact your dealer.

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